

TN-Advanced™

Tibial Nailing System

Surgical Technique

Suprapatellar Approach

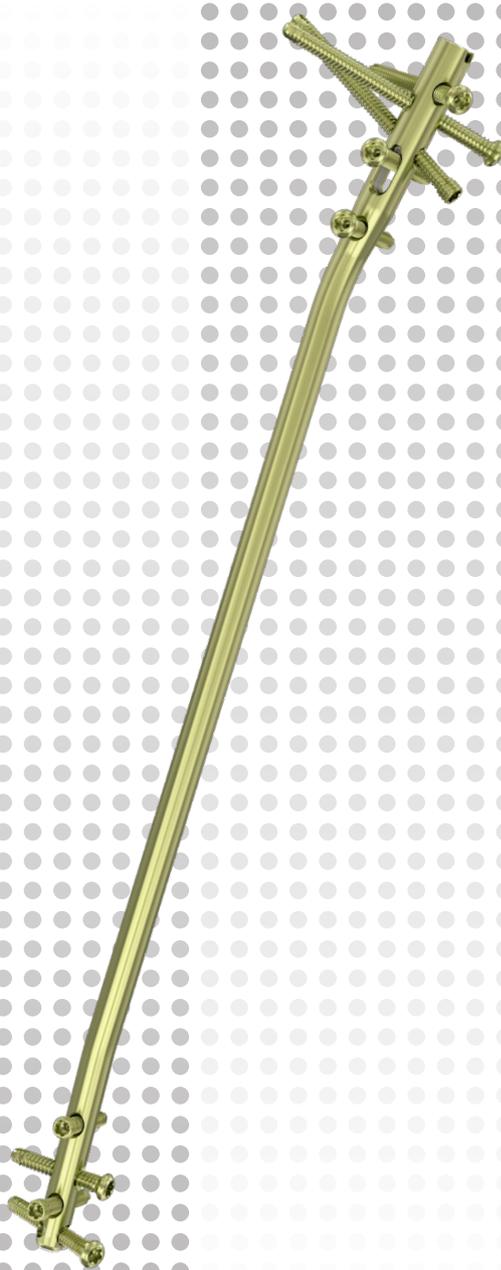


 Image intensifier control

This description alone does not provide sufficient background for direct use of DePuy Synthes products. Instruction by a surgeon experienced in handling these products is highly recommended.

Processing, Reprocessing, Care and Maintenance

For general guidelines, function control and dismantling of multi-part instruments, as well as processing guidelines for implants, please contact your local sales representative or refer to:

<http://emea.depuysynthes.com/hcp/reprocessing-care-maintenance>

For general information about reprocessing, care and maintenance of DePuy Synthes reusable devices, instrument trays and cases, as well as processing of DePuy Synthes non-sterile implants, please consult the Important Information leaflet (SE_023827) or refer to:

<http://emea.depuysynthes.com/hcp/reprocessing-care-maintenance>

Table of Contents

Introduction	The AO Principles of Fracture Management	2
	Warnings and Notes	3

Surgical Technique	Opening the tibia	4
	Nail insertion	17
	Distal locking	31
	Proximal locking	38
	End Cap Insertion	51
	Implant Removal	52

Product Information	Implants	59
	Alternative Instruments	61
	Instruments	63
	Optional Instruments - Suprapatellar	70
	Optional Instruments	71

The AO Principles of Fracture Management

Mission

The AO's mission is promoting excellence in patient care and outcomes in trauma and musculoskeletal disorders.


AO Principles^{1,2}

1.



Fracture reduction and fixation to restore anatomical relationships.

2.

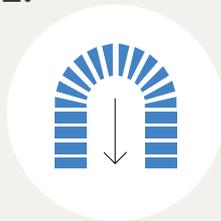

Fracture fixation providing absolute or relative stability, as required by the “personality” of the fracture, the patient, and the injury.

3.



Preservation of the blood supply to soft-tissues and bone by gentle reduction techniques and careful handling.

4.



Early and safe mobilization and rehabilitation of the injured part and the patient as a whole.

¹ Müller ME, Allgöwer M, Schneider R, Willenegger H. Manual of Internal Fixation. 3rd ed. Berlin, Heidelberg New York: Springer 1991.

² Buckley RE, Moran CG, Apivatthakakul T. AO Principles of Fracture Management: 3rd ed. Vol. 1: Principles, Vol. 2: Specific fractures. Thieme; 2017.

Warnings and Notes

▲ WARNINGS:

- Physician should take into account an increase in medullary pressure occurring during medullary nailing or reaming. This releases varying amounts of bone marrow and fat into the venous blood system.
- Physician should consider patient bone quality to ensure it provides adequate fixation to promote healing.

■ Notes:

- It is critical to ensure proper selection of the implant meets the needs of the patient anatomy and the presenting trauma. Physician should consider reaming to avoid under-sizing, to improve fit of nail, and to accelerate bone healing.
- Use of these devices is not recommended when there is systemic infection, infection localized to the site of the proposed implantation or when the patient has demonstrated allergy or foreign body sensitivity to any of the implant materials.
- Conditions that place excessive stresses on bone and implant such as severe obesity or degenerative diseases, should be considered. The decision whether to use these devices in patients with such conditions must be made by the physician taking into account the risks versus the benefits to the patients.
- Compromised vascularity in the site of proposed implantation may prevent adequate healing and thus preclude the use of this or any orthopaedic implant.

MRI Information on Torque, Displacement, Image Artifacts and Radio-Frequency-(RF-)induced heating can be found in the corresponding System Instructions for Use.

Please refer to the corresponding Instructions for Use for specific information on Intended use, Indications, Contraindications, Warnings and Precautions, Potential Adverse Events, Undesirable Side Effect and Residual Risks.

Instruction for Use are available at www.e-ifu.com and/or www.depuysynthes.com/ifu

Opening the tibia

1. Position patient

1A Position the patient supine on the radiolucent table. Ensure that the knee of the injured leg is flexed about 10° – 30° by using a knee roll or hydraulic leg holder. The knee roll can be placed under the lower part of the thigh if it obstructs the view of the tibial plateau in the AP view.

1B Position the image intensifier so that visualization of the tibia, including the articular surface proximally and distally, is possible in AP and lateral views.




2. Reduce fracture

2A Perform closed reduction manually by axial traction under image intensifier.

2B The use of a large distractor may be appropriate in certain circumstances.

■ **Note:**

The reduction can be temporarily fixed with reduction clamps. In epiphyseal fractures, the condyles or the pilon are fixed first in order to enable the nail insertion.

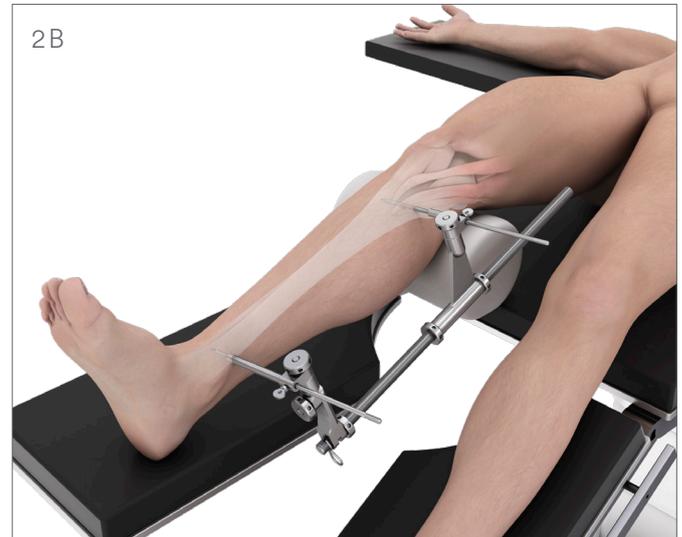

3. Make incision

3A Make a 2–4 cm longitudinal skin incision 4 cm proximal to the superior pole of the patella. The deep incision, also longitudinal, splits the quadriceps tendon in its midsubstance, just above its insertion into the patella and enters the knee joint through the suprapatellar pouch. Blunt dissection can be used to loosen the patella in the suprapatellar pouch, allowing the patella to lift off. Displace the patella anteriorly.

■ **Note:**

In arthritic knees, blunt release with an elevator or cauterization can be done in the suprapatellar pouch area medially and laterally. If these methods are not sufficient, extend the incision medially or laterally along the patella.



4. Determine entry point

Instruments

03.043.033S Protection Sleeve, for Nails, Ø 8–11 mm, flexible, long, sterile

or

03.043.035S Protection Sleeve, for Nails, Ø 8–11 mm, rigid, long, sterile

03.043.005 Wire guide, multihole, 8–11 mm and

03.043.008S Trocar, Ø 12 mm, long, sterile

or

03.043.034S Protection Sleeve, for Nails Ø 8–13 mm, flexible, long, sterile

or

03.043.036S Protection Sleeve, for Nails, Ø 8–13 mm, rigid, long, sterile

03.043.009 Wire guide, multihole, 8–13 mm and

03.043.012S Trocar, Ø 14 mm, long, sterile



03.045.018 Guide Wire Ø 3.2 mm, with drill tip, length 400 mm

4A In the frontal view the entry point is in line with the axis of the intramedullary canal and with the lateral tubercle of the intercondylar eminence.

4B In the sagittal view the entry point is at the ventral edge of the tibial plateau.

▲ Precaution:

Deviation from the optimal entry portal may cause irreducible malalignment, iatrogenic bone and soft tissue damage, malunion, and non-union.

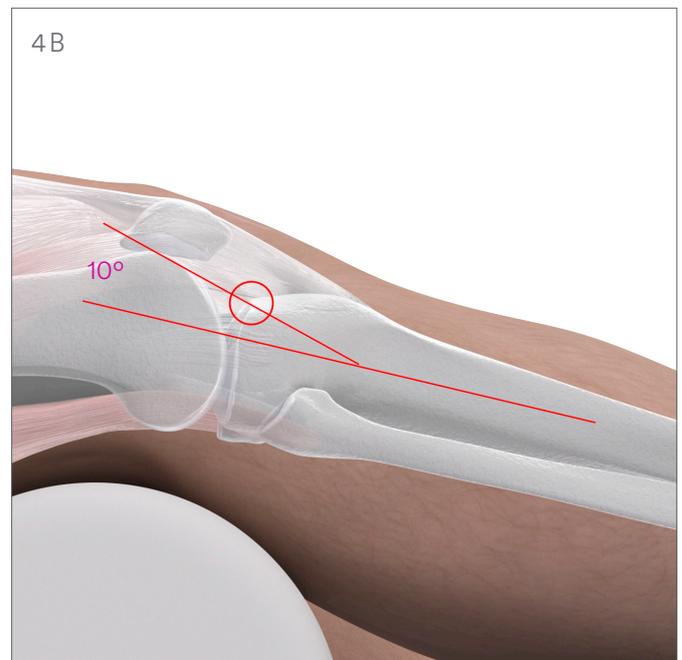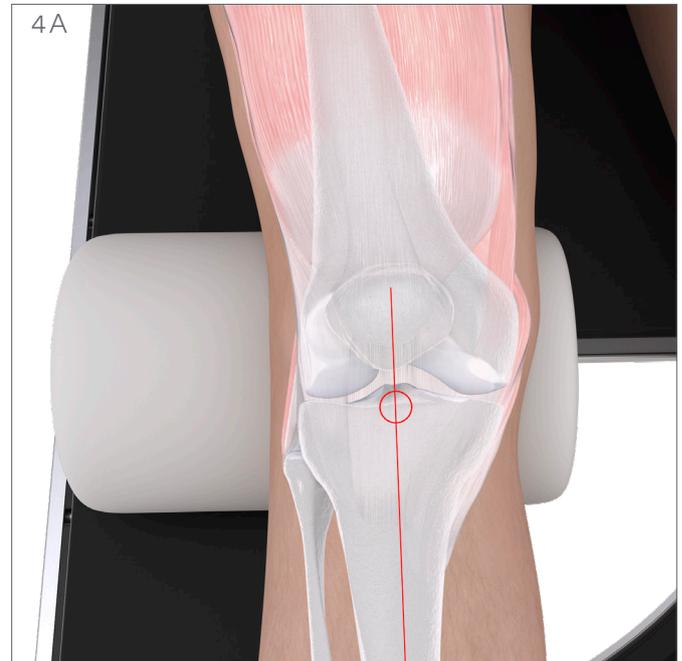


4C Assemble the trocar to the protection sleeve.


4D Insert the assembly through the incision into the knee joint, so that it glides between the articular surface of the patella and the trochlea of the distal femur and rests securely in this groove, while the patella is displaced anteriorly above the sleeve. Allow the trocar to slide out, as the protection sleeve is advanced to the anterior surface of the proximal tibia.



4E Unlock the trocar by rotating it a quarter turn and remove it from the sleeve.

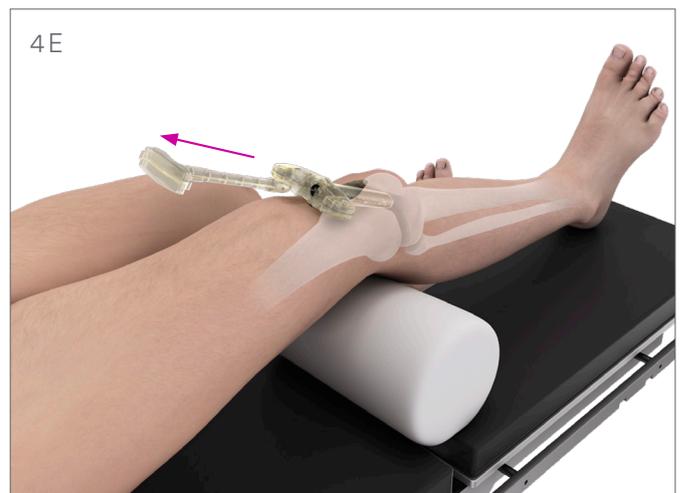

4F Insert the wire guide into the protection sleeve and advance it to the anterior surface of the proximal tibia.

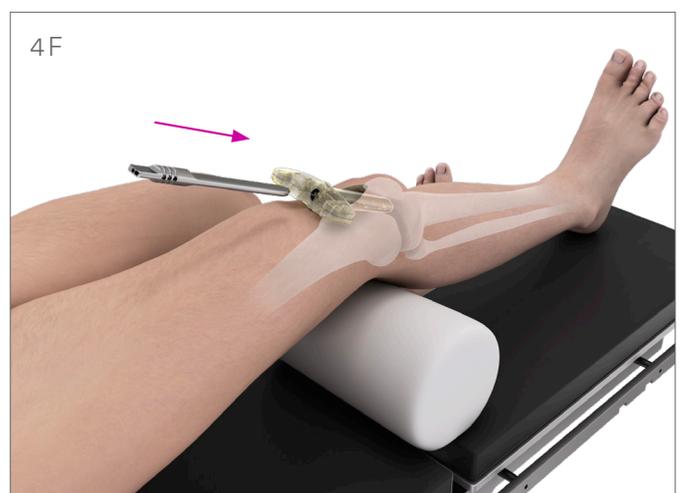


4G Insert a guide wire approximately 8 cm–10 cm and check the position under imaging in the AP and lateral views.

▲ Precaution:

Flexion of the knee must not be changed once the guide wire is inserted. A change could lead to increased pressure on the cartilage and could hinder surgical steps.



4H Adjustments to the guide wire location can be “dialed-in” by rotating the wire guide to place a second guide wire while the first guide wire remains in place.

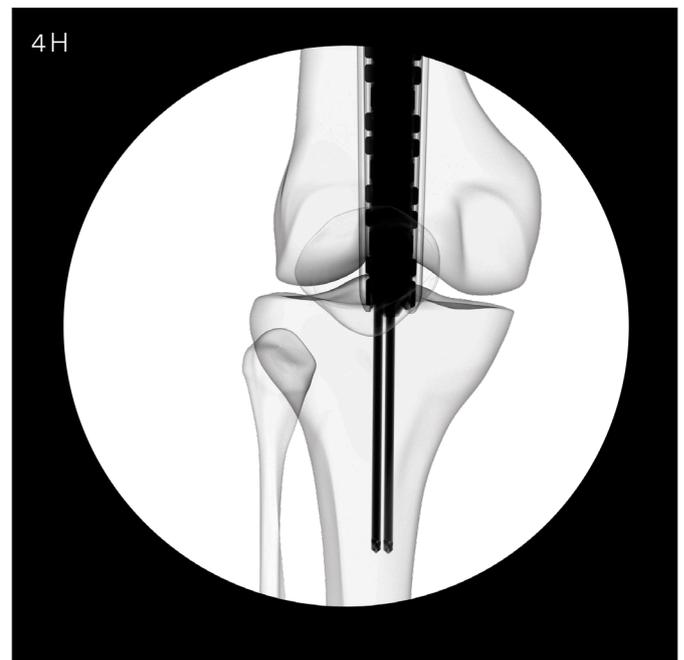

4I After correct placement of the second guide wire, remove the initial, central guide wire.



4J Take out the wire guide and re-insert it over the correctly placed guide wire through its central hole, to ensure that the instruments are centered over the correct entry point.



4. ALT.: Alternative Technique

Alternatively, the wire guide can be used to determine the entry point, before the insertion of the protection sleeve, thereby increasing maneuverability of the wire guide.

4. ALT. A Insert the wire guide through the incision into the knee joint, so that it glides between the articular surface of the patella and the trochlea of the distal femur and rests securely in this groove, while the patella is displaced anteriorly. Advance the wire guide to the anterior surface of the tibia. Slight adjustments of the knee flexion (between 10° and 30°) will provide the ideal radiographic location for the starting point and insertion of a guide wire.

4. ALT. B Insert the guide wire approximately 8 cm–10 cm and check the position under imaging in the AP and lateral views.

4. ALT. C Insert the protection sleeve over the wire guide and advance it to the anterior surface of the tibia.


▲ Precautions:

- Do not apply forces to the wire guide to reach the correct entry point. Adjust position of the wire guide by slight adjustments of the knee flexion (between 10° and 30°).
- Flexion of the knee must not be changed once the guide wire is inserted. A change could lead to increased pressure on the cartilage and could hinder surgical steps.


4. ALT. A



4. ALT. B

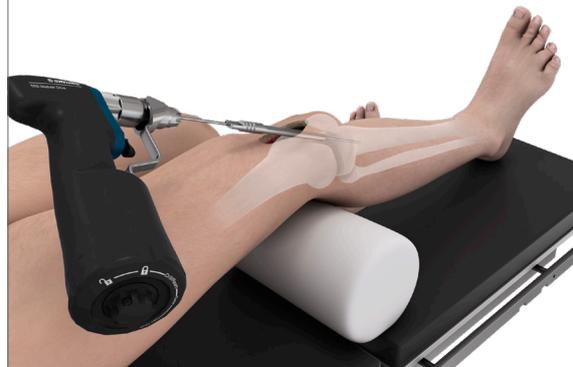


4. ALT. C

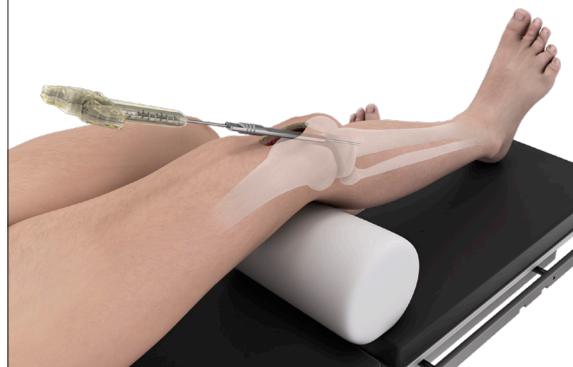

4. Option: Anchor the protection sleeve

4. OPT. A Once the protection sleeve is in its final position, it can be anchored to the femur using 3.2 mm Guide wire (e.g. 03.010.115).

4. OPT. A

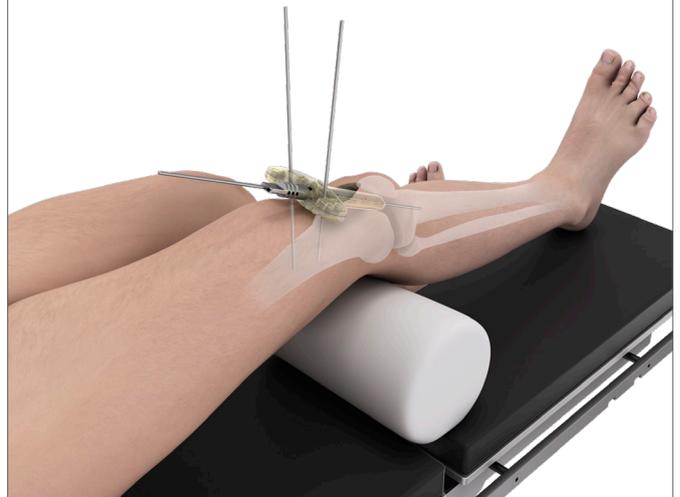

4. OPT. B If fixation to the tibial plateau is desired, the alternative sleeves 03.043.007S, 03.043.011S, 03.043.037S & 03.043.038S allow fixation to the tibial plateau, using guide wires 03.043.013.

4. OPT. B



5. Open medullary canal

Instrument

03.043.016 Drill Bit

5A Remove the wire guide from the protection sleeve and place the drill bit over the guide wire, through the protection sleeve, and down to the bone.



5B Drill to a depth of approximately 8–10 cm.

▲ Precaution:

Pay special attention not to penetrate the posterior cortex.

Remove the drill bit and guide wire. When using 12 and 13 mm nails, the opening must be over-reamed by at least 1 mm using a medullary reamer system.

▲ Precaution:


Do not drill inside the protection sleeve.

■ Note:

Dispose of the guide wire. Do not reuse.

▲ Precaution:

The suprapatellar protection sleeve is available in two different diameters. Markings on the sleeves indicate compatible nail diameters. Suprapatellar protection sleeves allow nail insertion through the sleeve and are compatible with SynReam reamer heads which are up to 1.5 mm larger in diameter than the largest compatible nail.

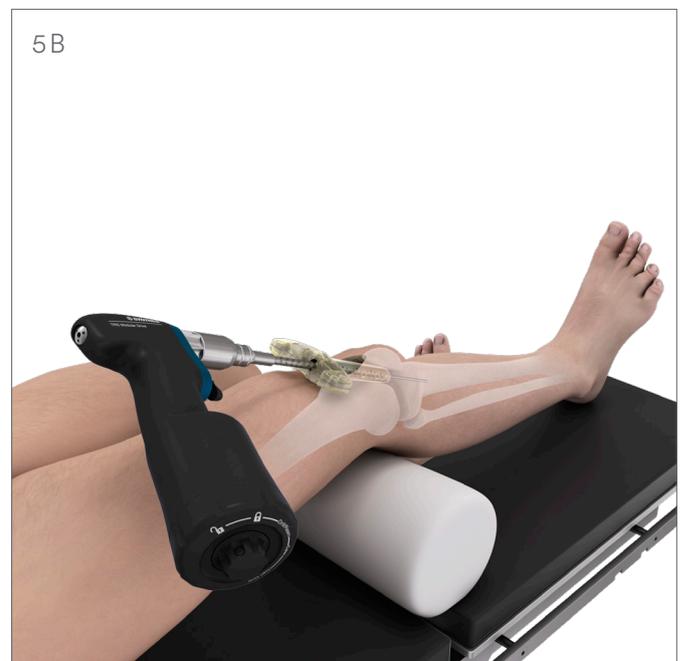

5. OPT.: Reaming Rod

Instruments

03.233.011S	Reaming Rod, Ø 3 mm, L 1150 mm or
351.704S	Reaming Rod, Ø 2.5mm, L 1150 mm, w/ball tip & extension
03.010.495	IMN Reduction Tool
03.010.496	T-Handle With Quick Coupling
03.010.093	Rod Pusher
03.045.035	Direct Measuring Device
03.045.036	Elongation Tube

5. OPT. A The use of a reaming rod can facilitate reduction, serve as a guide for intramedullary reamers, and aids in keeping bone fragments aligned during nail insertion.

▲ Precaution:

The Tibial Nail Advanced is cannulated and can be inserted over reaming rods with a diameter of up to 3.8 mm at their widest point. Compatible reaming rods will pass through the dedicated hole in the center of the aiming arm

5. OPT. B The IMN reduction tool may be used to aid in achieving alignment of the proximal and distal fragments, and for guiding the reaming rod into the distal fragment.

Insert the reduction finger to the desired depth.

5. OPT. A



5. OPT. B




5. OPT. C Pass the reaming rod through the cannulation of the instrument.

Remove the reduction instrument.

■ **Note:**

Use the rod pusher to help retain the reaming rod during the extraction of the reduction instrument.

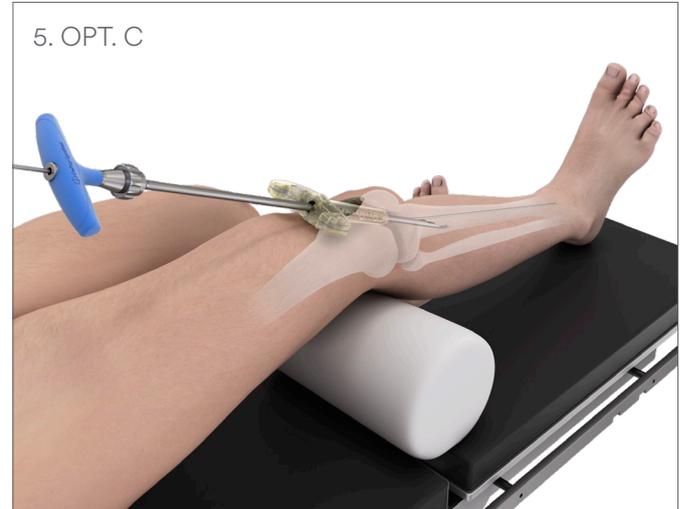

5. OPT. 2: Ream the medullary canal

Instruments

03.233.011S	Reaming Rod, Ø 3 mm, L 1150 mm or
351.704S	Reaming Rod, Ø 2.5mm, L 1150 mm, w/ball tip & extension

5. OPT. 2A If deemed appropriate, enlarge the medullary canal to the desired diameter using a DePuy Synthes reamer system intended for tibial reaming procedures by following the corresponding instructions of the reamer system.

- Protection sleeves for nails 8–11 mm allow reamer heads up to 12.5 mm.
- Protection sleeves for nails 8–13 mm allow reamer heads up to 14.5 mm.

5. OPT. 2A Check fracture reduction under image intensifier control.

▲ **Precautions:**

- Do not ream inside the protection sleeve.
- Monitor that the tip of the protection sleeve remains in direct contact to the proximal tibia.



▲ Precaution:

The reamer must travel through the protection sleeve before entering the bone. This may require a longer reamer shaft.

5. OPT. 2C We recommend using a long reamer shaft for nails exceeding 315 mm in length, in order to allow reaming along the entire length of the nail.

Option: Use the rod pusher (03.010.093) to prevent the reaming rod from backing off.

6. Determine Nail Length

Instruments

03.045.035	Direct Measuring Device
03.045.036	Elongation Tube

6A The nail length is determined either radiographically or over the reaming rod.

6B To measure nail length over the reaming rod, confirm reaming rod insertion depth under image intensification and account for a possible distraction or shortening at the fracture site.

6C Assemble the direct measuring device and elongation tube and pass the assembly over the reaming rod and down to the nail entry point.

6D Read the nail length directly from the measuring device.

6E The measuring scale is calibrated to the total length of a 950 mm reaming rod. When using a 1150 mm reaming rod, the length is indicated by the etched line on the reaming rod.

6F The nail diameter is determined either by reaming (optional) or radiographically.




Nail insertion

1. Assemble insertion instruments

Instruments

03.043.024	Insertion Handle
03.043.025	Connecting Screw
03.043.027	T-Handle Ball Hex Screwdriver

1A Connect the hexagonal ball at the tip of the screwdriver to the recess of the connecting screw by pushing both parts together until they snap into place and the connecting screw is retained to the screwdriver.

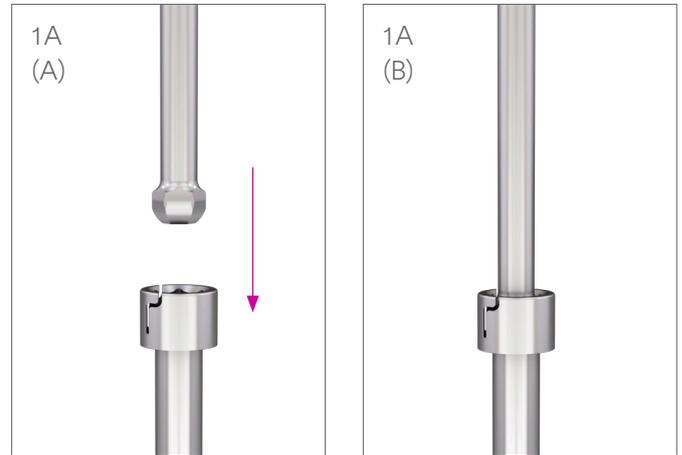

1B Connect the insertion handle to the nail by aligning the markings on the nail with the two slots on the barrel of the insertion handle. Push both parts together until they snap into place. The connection holds the nail in place until the connecting screw is tightened.





1C Pass the connecting screw through the insertion handle to engage with the nail and securely tighten it with the screwdriver. Remove the screwdriver.

▲ Precautions:

- Ensure that the connection between the nail and the insertion handle is tight. Retighten if necessary, after hammering and prior to the attachment of the aiming arm.
- Do not attach the aiming arm to the insertion handle at this point.

■ Note:

If using 03.043.035s, 03.043.036s, 03.043.037s or 03.043.038s. The inner metal sleeve must be removed prior to nail insertion.

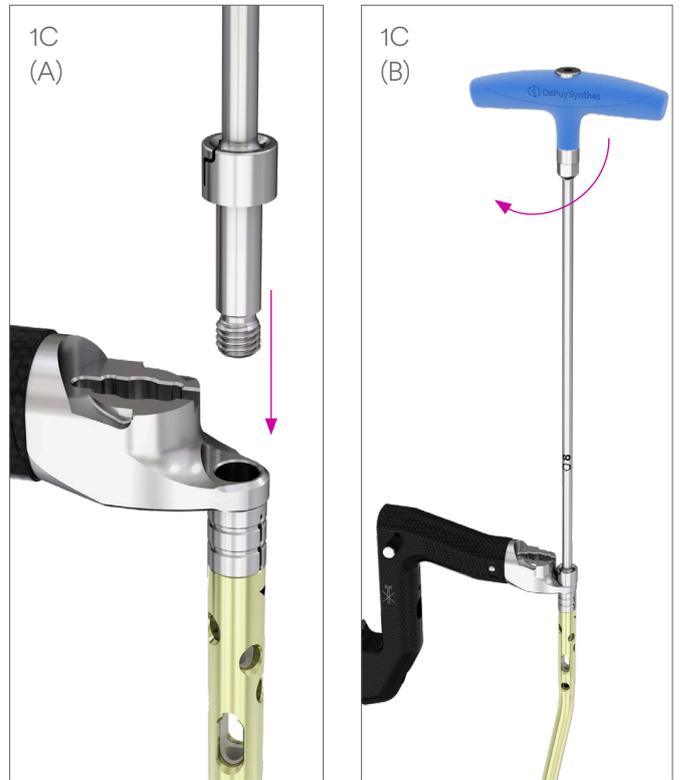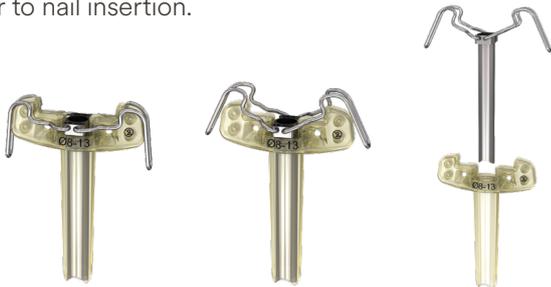

2. Insert Nail

Instruments

03.010.522	Hammer
03.043.028	Driving Cap
03.043.027	T-Handle Ball Hex Screwdriver
03.010.170	Hammer Guide

2A Insert the nail into the intramedullary canal through the protection sleeve.

2B Monitor the nail passage across the fracture; control in two planes to avoid malalignment.

2C Insert the nail until it is at or below the tibial opening. Check final nail position in AP and lateral views.



▲ Precaution:

To use the hammer, attach the driving cap to the insertion handle and secure it by twisting it one quarter turn. Apply light and controlled hammer blows to the driving cap to seat the nail.

■ Notes:

- The hammer guide may aid in controlling the direction of the hammer blows. Therefore, the hammer guide can be attached to the back end of the drive cap by screwing both parts together.
- After using the hammer, retighten the connecting screw to the nail with the screwdriver.
- If the nail needs to be backed out with a hammer, use the slotted hammer to slide along the bar of the driving cap or the hammer guide. Do not strike directly on the insertion handle as this will damage the handle.

▲ Precautions:

- If insertion is difficult, use the c-arm to confirm that there is no obstruction of the medullary canal. If no obstruction is found, choose a nail with a smaller diameter or enlarge the entry canal by reaming the medullary canal to a larger diameter.
- Do not use excessive twisting motions of the insertion handle.
- Remove Reaming Rod.



3. Check proximal nail position


Instruments

03.043.029	Aiming Arm
03.045.018	Guide Wire


3A To check the insertion depth of the nail, insert a 3.2 mm guide wire through the hole in the insertion handle as shown in the illustration.

3B Check proximal nail position under image intensifier control in the lateral view.

The tip of the guide wire indicates the exact proximal position of the tibial nail.

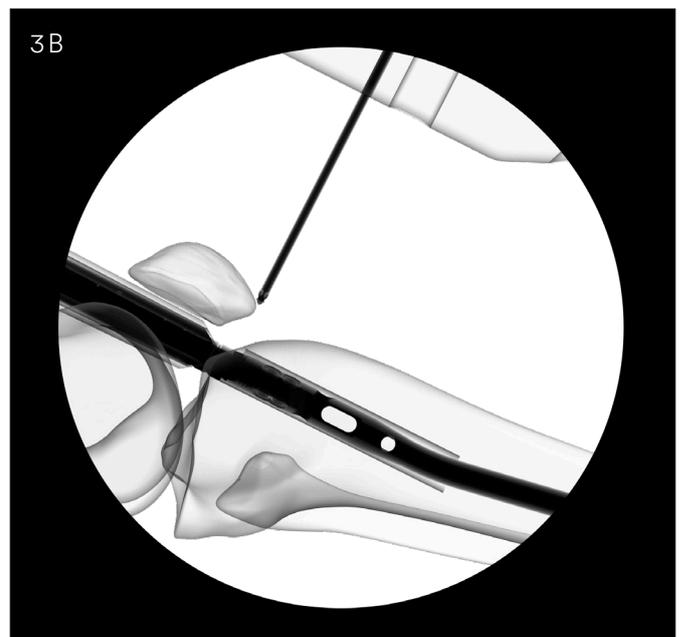


3C At this point, the trajectories of the three most proximal locking screws can be projected on a c-arm image. Attach the aiming arm to the insertion handle, by sliding it into the hook at the distal part of the insertion handle (1) and then rotating the latch towards the insertion handle for both parts to connect (2).

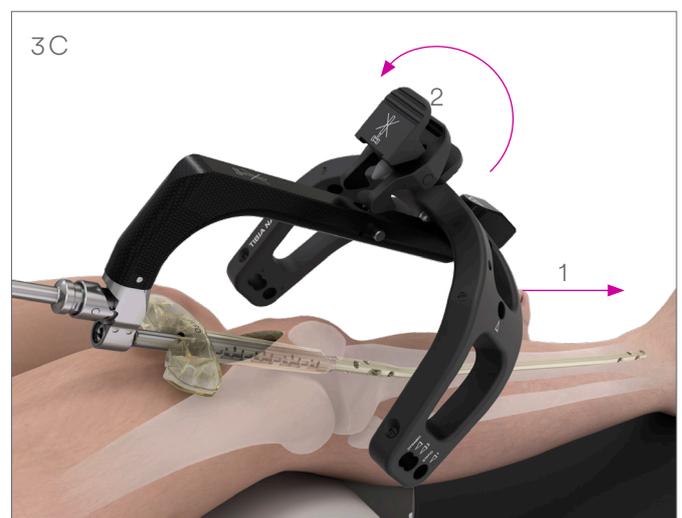

▲ Precautions:

- The distance between the markings on the insertion handle is 5 mm and corresponds to the extensions of the end caps. This feature can be used for over-insertion of the nail or for correcting the nail location within the medullary canal.
- If primary compression or secondary dynamization is planned, it is recommended to over-insert the nail by at least 7 mm, which corresponds to the maximum distance between the positions in static and dynamic modes. Protrusion of the proximal end of the nail can lead to irritation of the patellar tendon.

3D The trajectories of the three most proximal locking screws can be projected on a c-arm image by placing the drill bit 03.045.022 through the dedicated holes in the aiming arm. Insert the protection sleeve 03.045.019 and drill sleeve 03.045.020 into the corresponding hole in the aiming arm and assess the trajectory of the screw by taking a c-arm image in which the projections of the drill bit and the drill sleeve overlay.

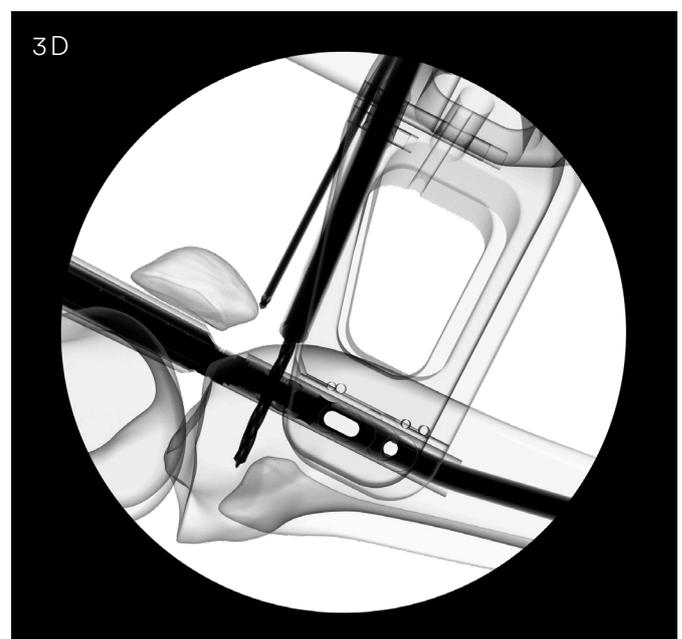

3E If the driving cap has been used, remove it now.

Remove the aiming arm, unless proximal locking is the next step.



4. Check distal nail position

Instrument

03.043.027 T-Handle Ball Hex Screwdriver

4A Check final nail position under image intensifier control in AP and lateral views.

4B Ensure that the reaming rod has been removed.


▲ Precaution:

Insertion depth is critical for distal third fractures where a minimum of two locking screws below the fracture line are required to stabilize the distal segment.

4C Confirm that the nail is securely connected to the insertion handle, especially after hammering, using the screwdriver.

▲ Precautions:

- To achieve compression the tibial nail needs to be locked distally first. The tibial nail allows a maximum compression or dynamization of 7 mm.
- Depending on the fracture patterns it might be advantageous to lock proximally first.

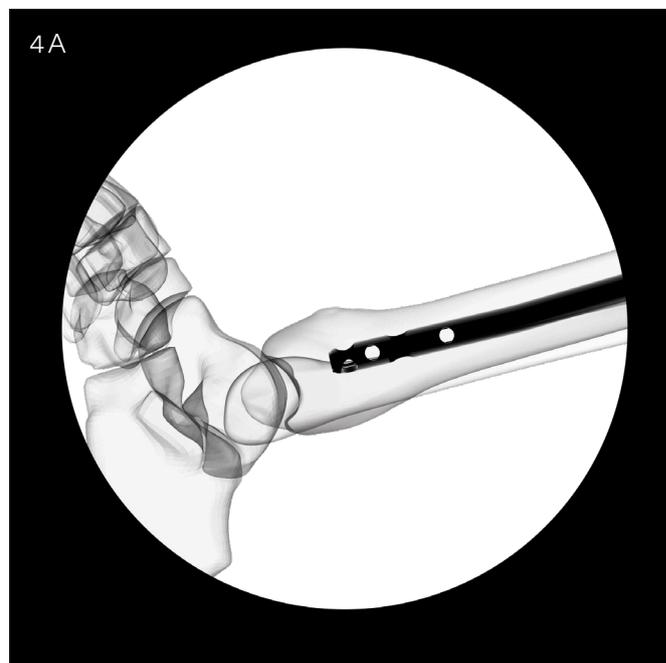


Locking

Locking Holes

Inlays made of polyetheretherketone (PEEK) in the nail provide angular stable locking in all screw holes, except for the Dynamic and Proximal Static 2 locking holes. The Ø 8 mm nail does not have a locking polymer inlay in the distal end of the nail.



Distal Oblique Locking Hole
Most distal screw offers choice of 30° lateral or medial angulation.

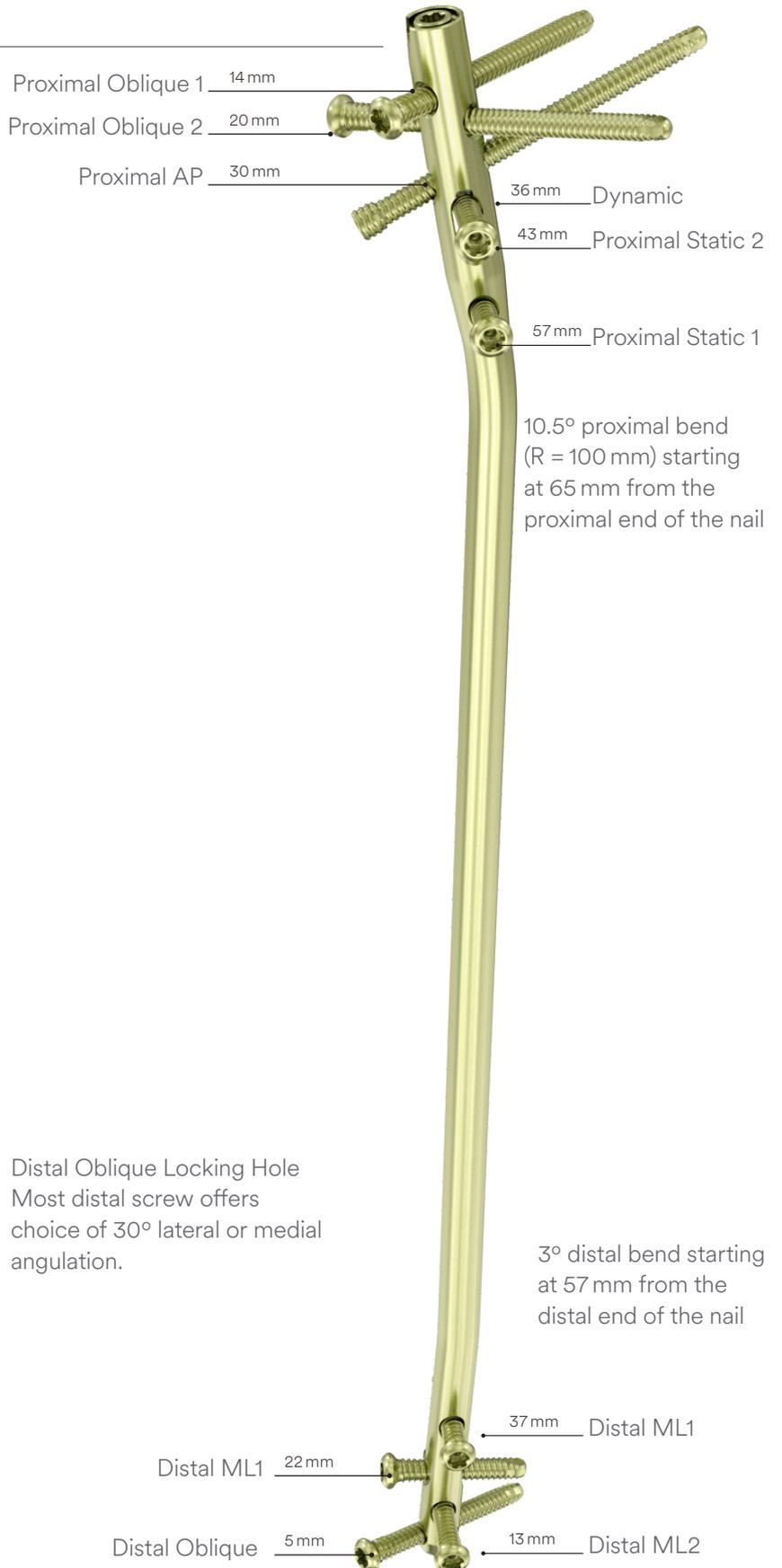

Locking


Use the appropriate locking screws and drill bit for the nail diameter selected.

Nail Diameter	Locking Screw	Drill Bit
9 mm to 13 mm	5.0 mm (light green)	4.2 mm
8 mm	Proximal Locking: 5.0 mm (light green)	4.2 mm
	Distal Locking: 4.0 mm (dark blue)	3.2 mm



Distal and proximal locking of the nail may be performed in either order.

If using a backstroke technique to reduce fracture gaps, distal locking must be performed prior to proximal locking. The hammer guide is attached to the driving cap and insertion handle. Light reverse hammer blows may be used to compress the fracture; monitor reduction radiographically.

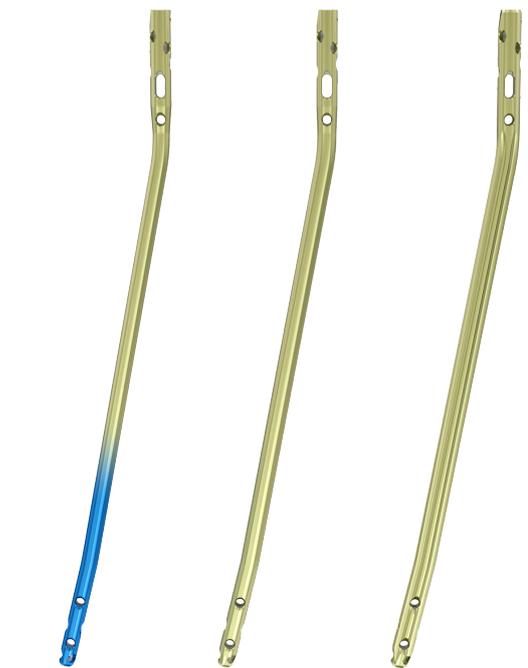

Ø 8 mm

Ø 9 mm
Ø 10 mm


Ø 11 mm
Ø 12 mm
Ø 13 mm

About Measuring Screw Length

Screw length is measured by using either of two methods.

1. Read length from the calibrated drill bits
2. Measure length using depth gauge for locking screws

Readings do not reflect the measured distance, they indicate the required screw length. The reading on the scale will correspond to the screw length as indicated on the screw label, taking into account the amount of screw tip protrusion required to get full screw thread engagement in the far cortex.

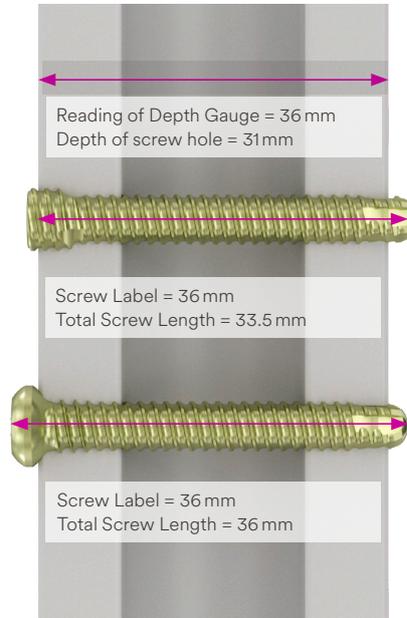

■ Notes:

- Drill bit location with respect to the far cortex is critical for measuring the appropriate locking screw length.
- Beware depth gauges are implant specific. Always use the appropriate depth gauge as specified in the surgical technique guide.

▲ Precaution:

Select adequate screw lengths to avoid protrusion of the screw tips and irritation of soft tissue.



Screw options

Tibial Nail Advanced offers two different types of screws:

1. Locking Screw
Standard IM nail locking screw.
2. Low Profile Locking Screw

■ Note:

- Both types of screws have a threaded recess and can be securely attached to the screwdriver by using the retention pins. To do so, slide the retention pin through the back of the screwdriver until it stops. Further advance it by turning it clock-wise, until its tip extends out of the tip of the screwdriver.
- The 4 mm screws do not have a threaded recess for attaching to the screwdriver with retention pins.

1. Locking Screw
2. Low Profile Locking Screw



Retention Pin



Engage the screwdriver in the recess of the locking screw and thread the retention pin into the screw's recess to lock the screw to the screwdriver.

4 mm Standard Screw



4 mm Low Profile Screw



Alternatively, the screw can be partially inserted with a power tool, by using the screwdriver shaft with its retention pin, following the same steps as described above.




▲ Precaution:

The screw must not be tightened with the power tool. Disengage the power tool from the screwdriver shaft before the screw is fully seated and use the manual handle to bring the screw to its final position.

■ Note:

Both types of screws have a threaded recess and can be securely attached to the screw-driver by using the retention pins (1). To do so, slide the retention pin through the back of the screwdriver until it stops (2.). Further advance it by pushing it through, until it clicks into place and its tip extends out of the tip of the screwdriver (3.).

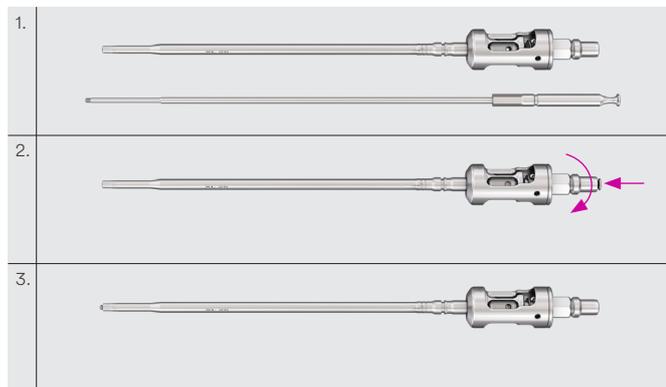

4. OPT.: Low Profile Screw

4. OPT. A The low profile screw can be used instead of the standard locking screw, by following the same basic steps for screw insertion.

4. OPT. B An optional sleeve is available to indicate when the screw is fully seated. Slide it over the tip of the screwdriver until it locks in place.


4. OPT. C In its initial position, it will cover the head of the screw, protecting surrounding soft tissues from the screw head's cutting flutes. Advance the screw until the sleeve touches the cortex.

■ Note:

Pay attention not to damage the cortex with the sleeve.

4. OPT. D Then retract the sleeve by pushing the release button and pulling it backwards towards the screwdriver handle.


4. OPT. A
Low Profile
Locking Screw

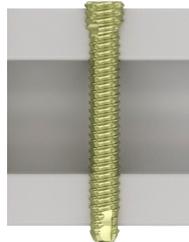


4. OPT. B



4. OPT. C

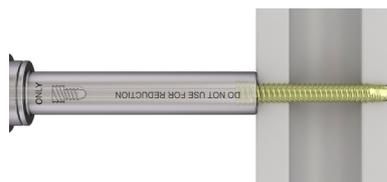

4. OPT. D

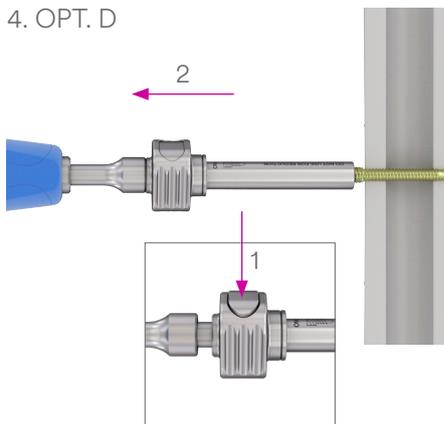

Distal locking

For distal third fractures a minimum of two locking screws are required to stabilize the distal segment, and the use of the most distal locking option is recommended. This locking option is oriented 30° from the sagittal plane.




The most distal screw hole allows placing the screw at an angle of 30° from the sagittal plane in either medial or lateral direction.




1. Drill for distal locking

Nail Diameter	Locking Screw	Drill Bit
9 mm to 13 mm	5.0 mm (light green)	4.2 mm 03.010.104
8 mm	Distal Locking: 4.0 mm (dark blue)	3.2 mm 03.010.103

- 1A Align image: Check the reduction, correct alignment of the fragments, and leg length before locking the nail. Align the c-arm with the hole in the nail closest to the fracture until a perfect circle is visible in the center of the screen. (Distal ML hole shown in illustration).

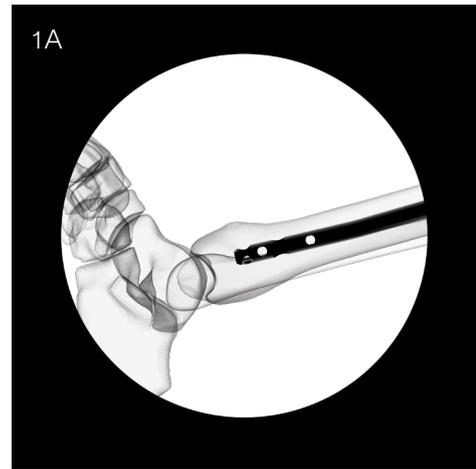


1B Determine incision point: Place a guide wire on the skin over the center of the hole to mark the incision point and make a stab incision. Insert the tip of the appropriate drill bit through the incision and down to the bone.

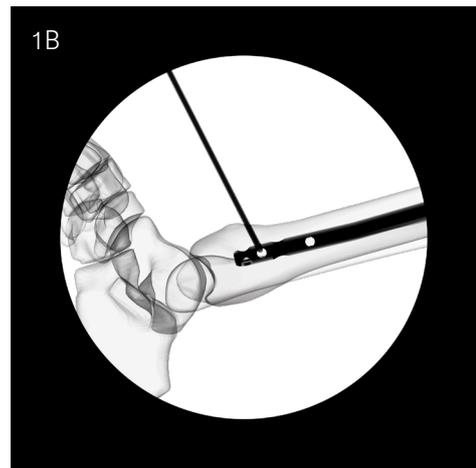

- 1C Drill: Use the image intensifier to position the tip of the drill bit over the center of the locking hole. Incline the drill bit so that it aligns with the axis of the beam of the image intensifier. Hold the drill bit in this position and drill through both cortices.

▲ Precaution:

Stop drilling immediately after the drill bit passes through the second cortex.

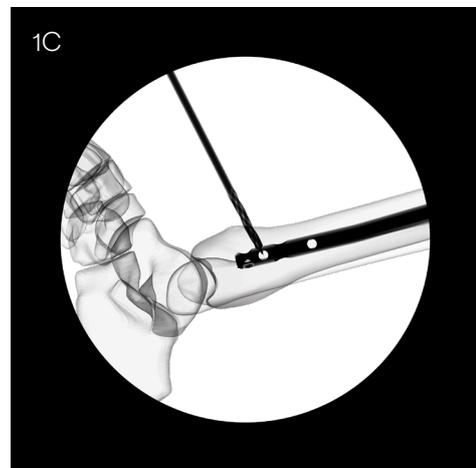


1. ALT.: Alternative instrument to drill for distal locking

Instrument

511.300 Radiolucent Drive (separately available)

Nail Diameter	Locking Screw	Drill Bit
9 mm to 13 mm	5.0 mm (light green)	4.2 mm 03.010.101
8 mm	Distal Locking: 4.0 mm (dark blue)	3.2 mm 03.010.100

1. ALT. A Alternatively, a radiolucent drive may be used for distal locking, following the same basic steps. The radiolucent drive enables the use of the image intensifier to ensure proper alignment of the drill with the locking hole in the nail. The radiolucent drive requires a separate set of drill bits.

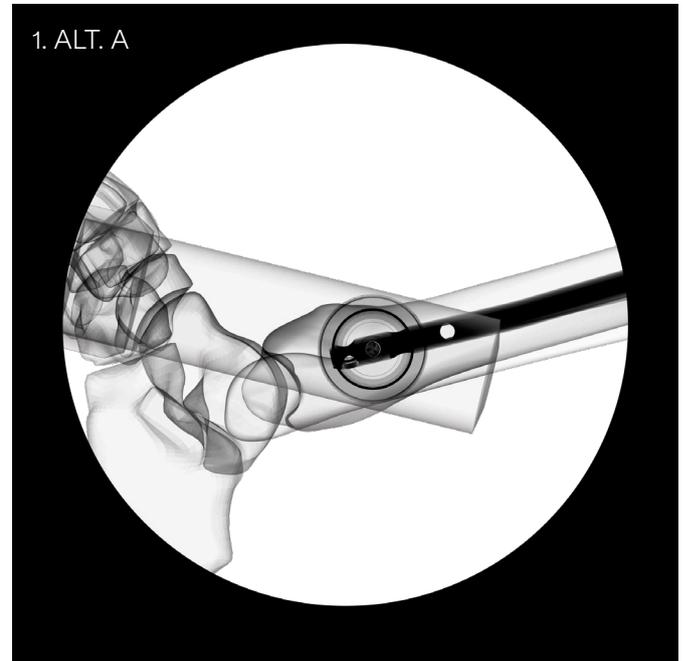


2. Determine locking screw length for distal locking

Instrument

03.010.429 Direct Measuring Device

2A Stop drilling immediately after both cortices and detach the drill bit from its coupling.

2B Under image intensifier control, ensure the correct position of the drill bit beyond the far cortex.

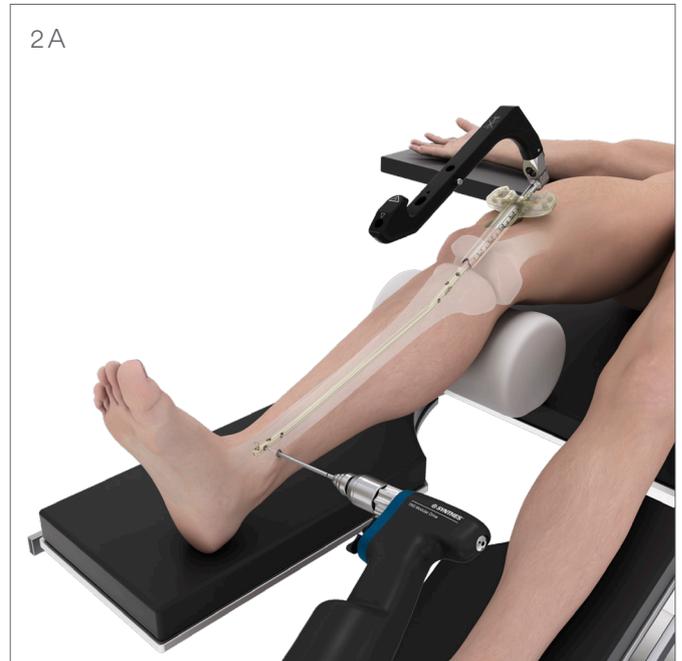

2C Place the direct measuring device onto the drill bit. Read the graduation of the measuring device at the end of the drill bit. This corresponds to the appropriate locking screw length.



2. ALT.: Alternative instrument to determine locking screw length for distal locking

Instrument

03.019.017 Depth Gauge

2. ALT. A Measure the locking screw length using the depth gauge for locking screws. Ensure the outer sleeve is in contact with the bone and the hook grasps the far cortex.

2. ALT. B Read the locking screw length directly from the depth gauge at the back of the outer sleeve.




3. Insert locking screw

Instruments

03.045.003 Screwdriver Short XL25

03.045.004 Retention Pin XL25

3A Insert the appropriate length locking screw using the screwdriver.

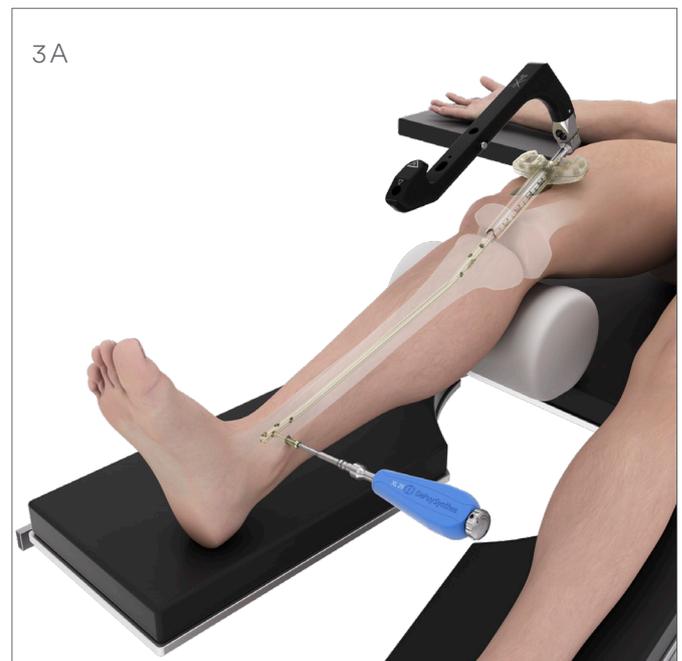

3. OPT.: Alternative instrument to insert locking screw

Instruments

03.045.007 Screwdriver Shaft Short XL25

03.045.008 Retention Pin


03.140.027 Handle With Quick Coupling

3. OPT. A Alternatively, the screw can be partially inserted with a power tool, by using the screwdriver shaft with its retention pin, following the same steps as described above.

▲ Precaution:

The screw must not be tightened with the power tool. Disengage the power tool from the screwdriver shaft before the screw is fully seated and use the manual handle to bring the screw to its final position.

3. OPT. A

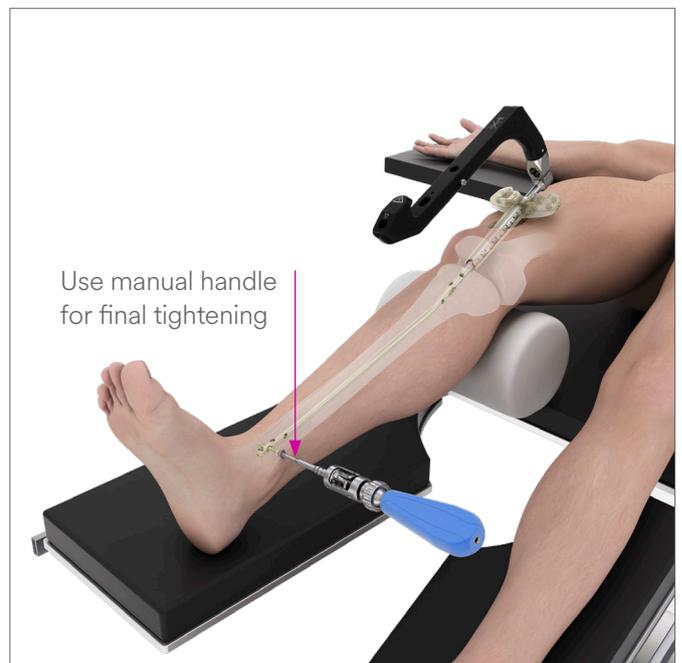

3. OPT. 2: Low profile locking screw

Instruments

03.045.010	Sleeve
03.045.029	Reamer Ø 5.5 mm

3. OPT. 2A The low profile screw can be used instead of the standard locking screw, by following the same basic steps for screw insertion.

3. OPT. 2B An optional sleeve is available to indicate when the screw is fully seated. Slide it over the tip of the screwdriver until it locks in place.

3. OPT. 2C The use of the Ø 5.5 mm reamer, to make room for the threaded screw head, is recommended in hard bone.

■ Note:


The use of the Ø 5.5 mm reamer, to make room for the threaded screw head, is a mandatory step for 4.0 mm low profile screws.

▲ Precaution:

Use a Ø 5.5 mm reamer, to make room for the threaded screw head of the 4.0 mm low profile locking screw for 8 mm nail.

3. OPT. 2A

Low Profile Locking Screw

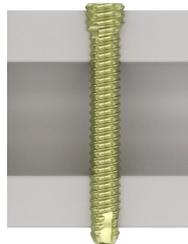


Proximal locking

Nail Diameter	Locking Screw	Drill Bit
All nails (8 mm to 13 mm)	5.0 mm (light green)	4.2 mm 03.045.022

▲ Precaution

(oblique and AP locking options): Proximal locking requires special attention. To avoid lesion of the popliteal artery, the tibial nerve and the common peroneal nerve, as well as damage to the proximal tibiofibular joint, drilling must be stopped immediately before penetrating the far cortex. Monitor the position of the drill bit.

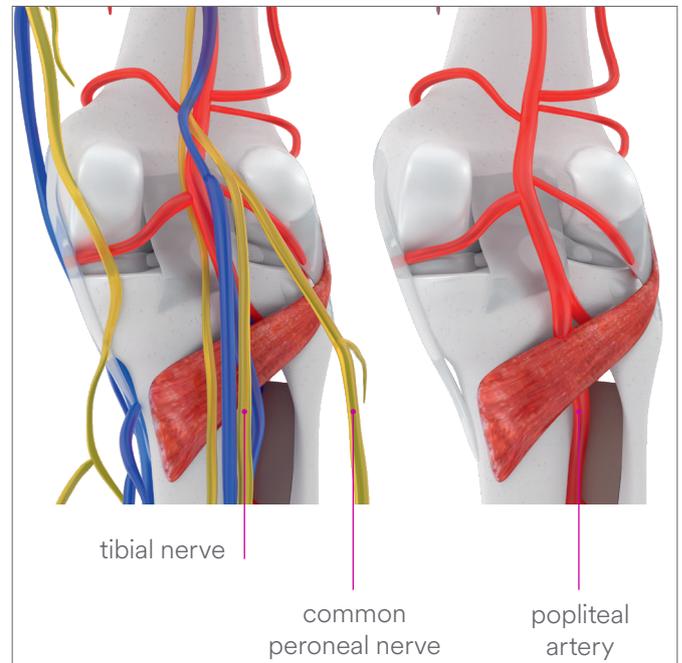

1. Mount aiming arm

Instruments

03.043.029	Aiming Arm
03.043.027	T-Handle Ball Hex Screwdriver

1A Confirm that the nail is securely connected to the insertion handle, especially after hammering, using the screwdriver.


1B Mount the aiming arm to the insertion handle.

▲ Precaution:

Do not exert forces on the aiming arm, protection sleeve, drill sleeves and drill bits. These forces may prevent accurate targeting through the proximal locking holes and damage the drill bits.

■ Note:

The proximal AP screw is inserted through the guiding hole in the insertion handle and does not require aiming arm to be attached.

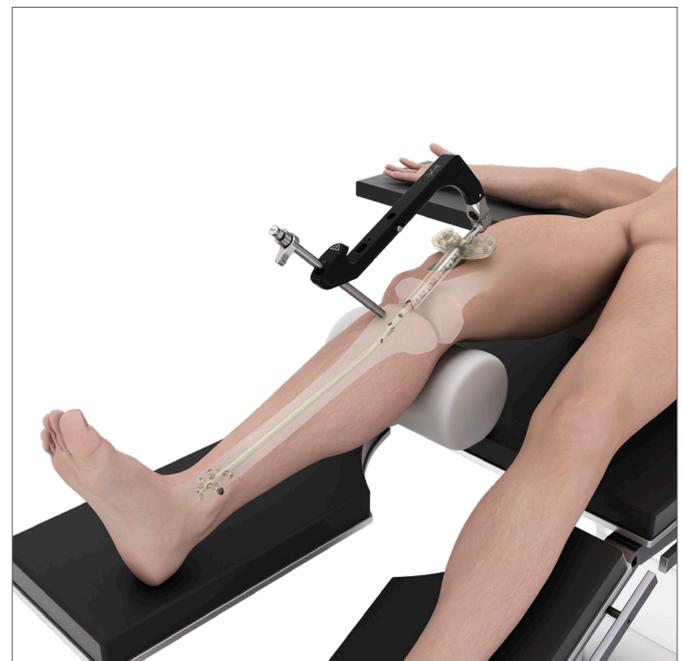

2. Insert trocar assembly

Instruments

03.045.019	Protection Sleeve
03.045.020	Drill Sleeve
03.010.070	Trocar Ø 4.2 mm

2A Assemble the three-part trocar assembly (protection sleeve, drill sleeve, and trocar). Align the triangular marking at the tip of the protection sleeve with the marking beside the desired hole on the aiming arm, and insert the three-part trocar assembly through the aiming arm. Make a stab incision and insert the trocar to the bone. Twist the protection sleeve by a quarter turn to lock it into place. Remove the trocar.




3. Drill for proximal locking

Instrument

03.045.022 4.2 mm Drill Bit

3A Ensure that the drill sleeve contacts the near cortex and insert the calibrated drill bit and start drilling the near cortex.

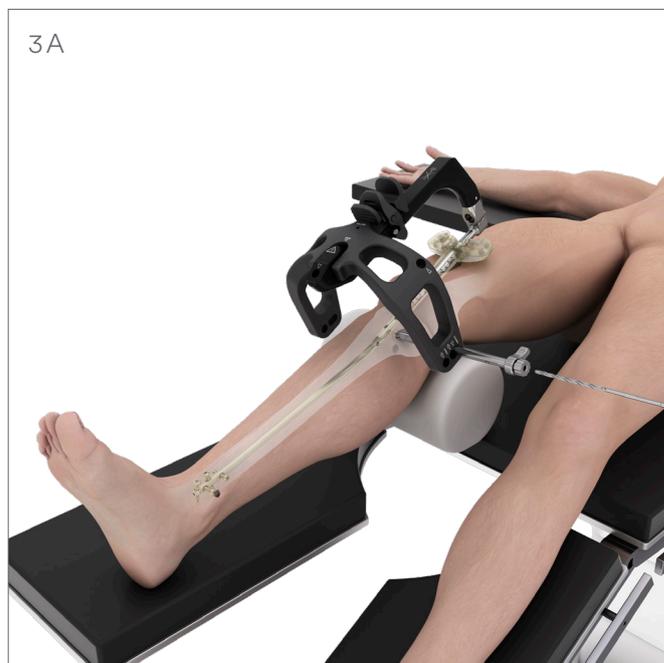

3B Drill to the desired depth and confirm drill bit position after drilling.

▲ Precaution

(medial to lateral locking options): Stop drilling immediately after penetrating both cortices.

▲ Precaution

(oblique and AP locking options): Proximal locking requires special attention. To avoid lesion of the popliteal artery, the tibial nerve, and the common peroneal nerve, as well as damage to the proximal tibiofibular joint, drilling must be stopped immediately before penetrating the far cortex. Monitor the position of the drill bit.



4. Measure screw length

4A Ensure that the drill sleeve contacts the bone and read the measurement from the calibrated drill bit at the back of the drill sleeve. This measurement indicates the appropriate length of the locking screw.



4B Remove the drill bit and the drill sleeve.




4. ALT.: Alternative technique

Instrument

03.019.017 Depth Gauge

4. ALT. A After drilling, remove the drill bit and the drill sleeve.

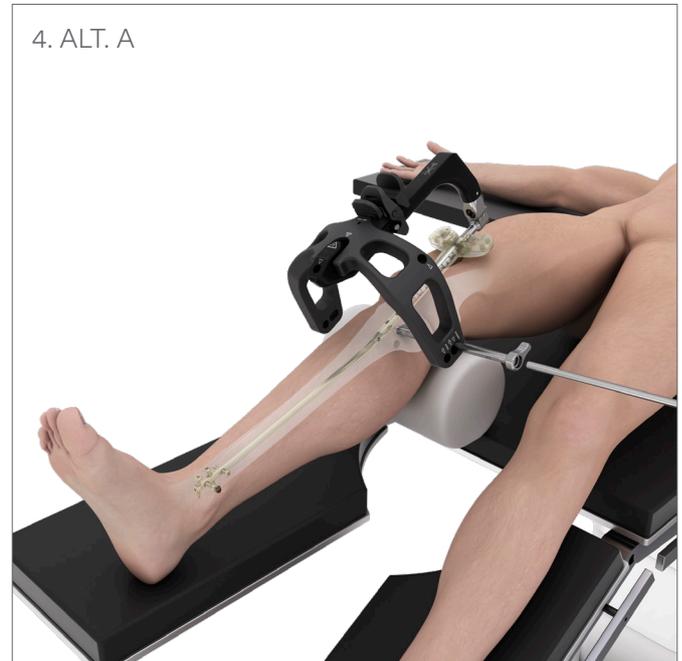

4. ALT. B Insert the depth gauge into the protection sleeve. Make sure that the hook grasps the far cortex for bi-cortical screws or touches the end of the screw hole for mono-cortical screws, and that its sleeve is on the bone.

Read the measurement from the back of the protection sleeve, which indicates the appropriate length locking screw.



5. Insert locking screw

Instruments

03.045.001 Screwdriver Long XL25

03.045.002 Retention Pin XL25

5A Insert the appropriate length locking screw using the screwdriver.



5. ALT.: Alternative instrument to insert locking screw

Instruments

03.045.005 Screwdriver Shaft Long XL25

03.045.006 Retention Pin

03.140.027 Handle With Quick Coupling

5. ALT. A Alternatively, the screw can be partially inserted with a power tool, by using the screwdriver shaft with its retention pin, following the same steps as described above.



▲ Precaution:

The screw must not be tightened with the power tool. Disengage the power tool from the screwdriver shaft before the screw is fully seated and use the manual handle to bring the screw to its final position.

5. ALT. A



5. ALT. B The shaft of the screwdriver has two lines, one of which indicates insertion depth of the locking screw (1), and the other indicating insertion depth of the low profile locking screw (2) relative to the tip of the protection sleeve. Screws are fully seated, when the line is flush with the head of the protection sleeve.



5. OPT.: Low profile locking screw

Instruments

03.045.009	Sleeve
03.045.029	Reamer Ø 5.5 mm

5. OPT. A The low profile screw can be used instead of the standard locking screw, by following the same basic steps for screw insertion.

5. OPT. B An optional sleeve is available to indicate when the screw is fully seated. Slide it over the tip of the screwdriver until it locks in place.

Before using the sleeve, unlock the protection sleeve, ensure it contacts the bone, and lock it into place again by twisting it a quarter turn.

5. OPT. C The use of the Ø 5.5 mm reamer, to make room for the screw head, is recommended in hard bone.

5. OPT. A



6. Option: Compression locking mode

Instruments


03.043.026	Compression Screw
03.043.027	T-Handle Ball Hex Screwdriver

6. OPT. A If the fracture gap needs compression after nail insertion, it can be accomplished without removing the insertion instruments. The nail allows a maximum compression of 7 mm.

6. OPT. B Distal locking is required prior to compression locking. Insert one proximal locking screw in the dynamic locking hole (DYNAMIC).

6. OPT. C Confirm that the nail is securely connected to the insertion handle using the screwdriver.

6. OPT. D Insert the compression screw through the connecting screw and into the nail using the screwdriver.


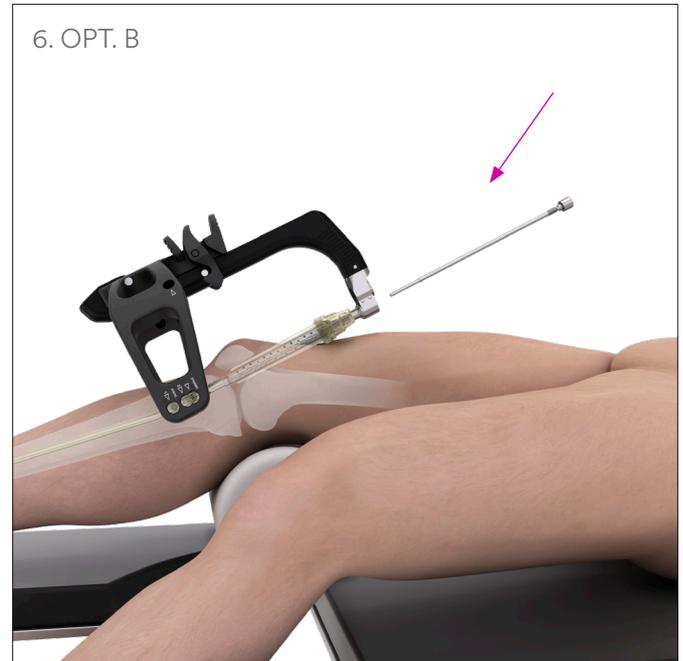


6. OPT. E The compression screw will contact the dynamic locking screw. Advance the compression screw until the fracture gap is reduced. Monitor reduction with the image intensifier.

6. OPT. F Each revolution of the compression screw corresponds to compression of 1 mm (maximum 7 mm). Control the fracture gap before, during, and after the compression procedure.

Note:

Do not overtighten the compression screw, it may deform the locking screw.

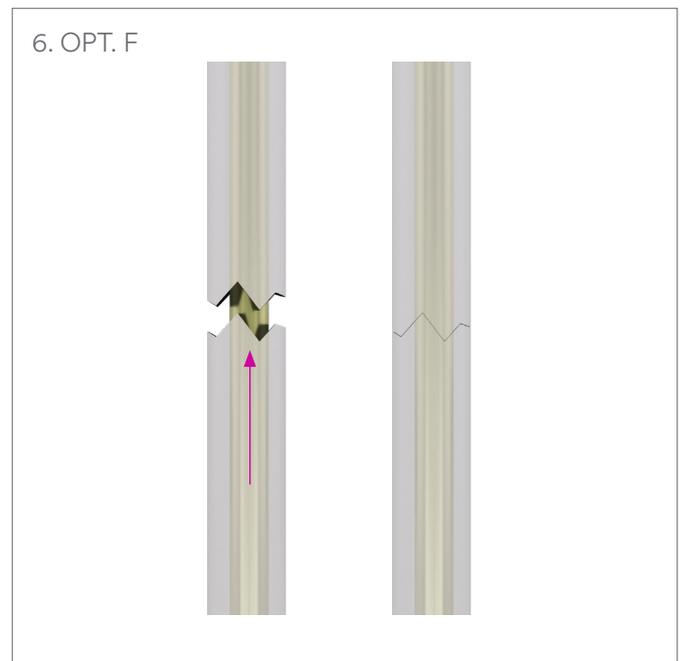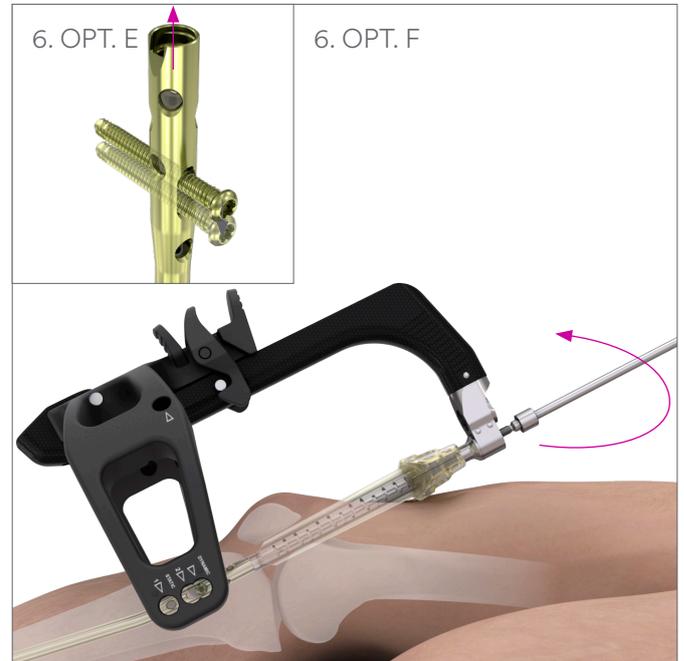

6. OPT. G Secure the compression by inserting a second proximal locking screw in the most distal hole of the proximal locking options (STATIC 1).

6. OPT. H Remove the compression screw.

6. OPT. I Additional oblique locking screw can be inserted if required.

6. OPT. G



6. OPT. H



End Cap Insertion

1. Insert the end cap

Instruments

03.045.001	Screwdriver Long XL25
03.045.002	Retention Pin XL25
03.043.027	T-Handle Ball Hex Screwdriver



1A Remove the connecting screw. The insertion handle can remain in place to help align the end cap to the top of the nail. The end cap fits through the barrel of the insertion handle.

▲ Note:

Flexion of the knee might have changed during the procedure. This can block the connecting screw and impede its removal. Slight adjustment of the knee's flexion (between 10° and 30°) will provide a neutral position that facilitates removal of the connecting screw and insertion of the end cap.

1B If desired, the end cap can be locked to the screwdriver by use of the retention pin. To do so, slide the retention pin through the back of the screwdriver until it stops. Further advance it by turning it clock-wise, until its tip extends out of the tip of the screwdriver.

1C Insert the end cap through the barrel of the insertion handle and tighten it to the nail.

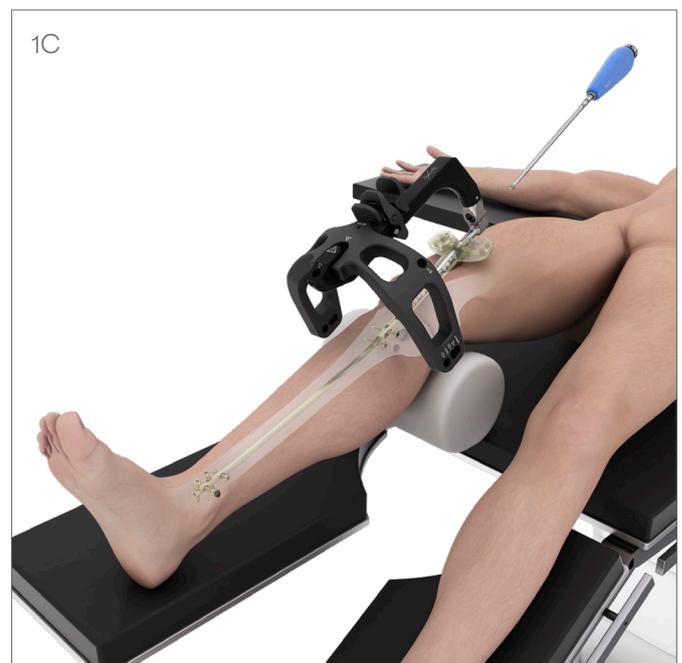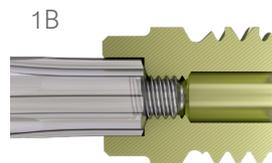

Implant Removal

Implant removal is an optional procedure.

The implant may be removed by a suprapatellar, or by a traditional infrapatellar approach.

If the implant is removed by a suprapatellar approach, the access portal needs to be prepared by following the same steps as described under steps 3 and 4. This includes the use of a suprapatellar protection sleeve.




If the implant is removed by an infrapatellar approach, position the patient supine on a radiolucent table. Ensure that the knee of the affected leg can be flexed at least 90°. Position the image intensifier such that visualization of the tibia proximally and distally is possible in AP and lateral views.

Make an incision in line with the central axis of the intra-medullary canal. Depending on the anatomy of the patient, this incision can be transpatellar, medial or even lateral parapatellar.

The incision starts proximally at the distal third of the patella along the patellar ligament down to the tibial tuberosity.

Mobilise the infrapatellar fat pad laterally and dorsally without opening the synovium, allowing access to the entry point of the nail.

Once the access portal has been made, the implant is removed by performing the following steps.

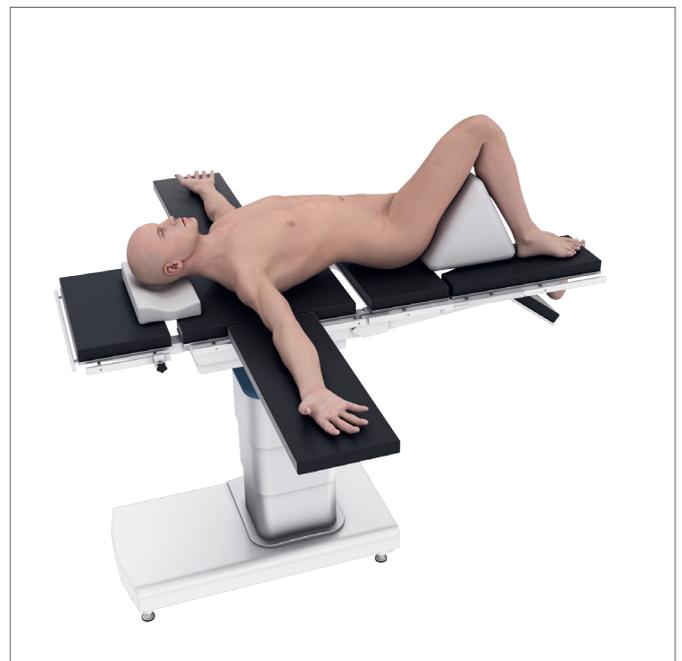

Remove end cap and locking screws

Instruments

03.045.001	Screwdriver Long XL25
03.045.002	Retention Pin XL25

Clear the recess of the end cap and the locking screws of any tissue ingrowth. Remove the end cap with the screwdriver.

■ Note:

The screwdriver (without the retention pin) will fit over a 1.6 mm guide wire (e.g. 292.655), which may be used to guide the screwdriver into the recess of the end cap. Once the tip of the screwdriver engaged with the recess of the end cap, the guide wire is removed through the back end of the screwdriver, and the retention pit can be used to secure the end cap to the screwdriver.

Remove all locking screws except one of the proximal locking screws with the screwdriver. Always remove the two most proximal locking screws in order to insert the extraction screw into the proximal end of the nail.

■ Note:

The XL25 recess is compatible with StarDrive SD25 screwdrivers.





Optional: Additional Instruments for Screw Removal

Instruments

03.045.030	Extractor Shaft for XL25 and SD25
03.045.031	Curette for XL25
03.045.032	Extraction Screw, conical
03.900.001	Sharp Hook, straight, length 150 mm

If screw heads are overgrown or the recess is damaged, additional instruments are available for screw removal. They can be used with all XL25 screw types.

1. Clear recess and screw head with the curette. The curette turns counter-clockwise
2. Use the sharp hook to clean out any remaining tissue
3. Engage the extractor shaft to remove the screw
4. If 3. does not work, use the conical extraction screw to remove the screw. The conical extraction screw turns counter-clockwise

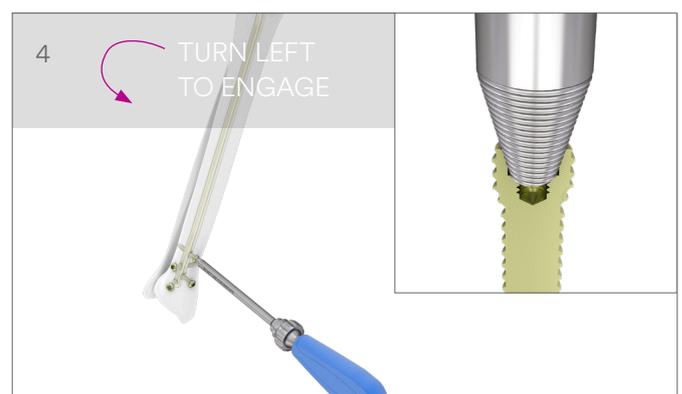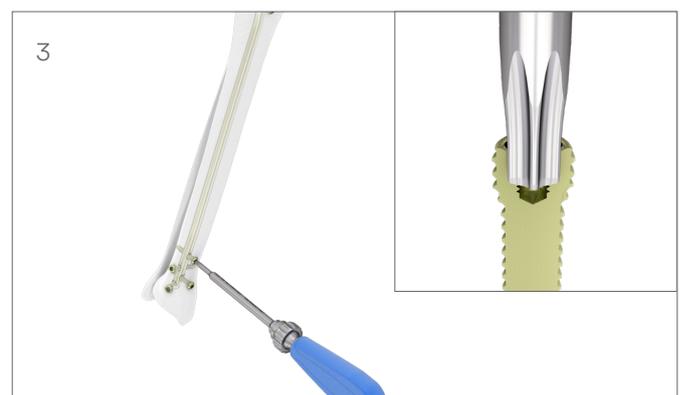

* Availability might be different depending on the market.

Attach extraction screw and hammer guide

Instruments

03.010.000	Extraction Screw, for Tibial and Femoral Nails
03.010.170	Hammer Guide
03.045.001	Screwdriver Long XL25
03.045.002	Retention Pin XL25

Before removing the final locking screw, screw the extraction screw into the tibial nail and tighten it to prevent rotation or displacement of the nail posteriorly below the tibial plateau.

Attach the hammer guide to the extraction screw.

Remove the remaining locking screw with the screwdriver.



Remove nail

Instrument

03.010.522 Hammer

Extract the nail by applying gentle blows with the hammer.



Alternative solid protection sleeve

Instruments

03.010.430	Handle for Protection Sleeve
03.010.435	Protection Sleeve 12.0, for 8.0–11.0 mm Nails
03.010.436	Protection Sleeve 14.5, for 12.0–13.0 mm Nails
03.010.437S	Outer Protection Sleeve 12.0, for 8.0–11.0 mm Nails, sterile
03.010.438S	Outer Protection Sleeve 14.5, for 12.0–13.0 mm Nails, sterile
03.010.455	Trocar Ø 12.0 mm
03.010.456	Trocar Ø 14.5 mm
03.010.433	Centering Sleeve 12.0 / 3.2
03.010.434	Centering Sleeve 14.5 / 3.2

Assemble the handle, inner protection sleeve, outer protection sleeve, and trocar. Rotate the knob to lock the assembly into the handle during insertion.

Insert the handle assembly through the incision into the knee joint, so that it glides between the articular surface of the patella and the trochlea of the distal femur and rests securely in this groove, while the patella is displaced anteriorly above the cannula. When the trocar reaches the surface of the tibia, rotate the knob on the handle to allow the trocar to slide out, as the cannula is advanced to the anterior surface of the proximal tibia. Remove the trocar. After removal of the trocar, insert the centering sleeve through the inner protection sleeve. Advance to the anterior surface of the tibia.

Continue at step **4G** of “Opening the tibia”, until “Nail insertion”.

■ Note:

Before nail insertion, the metal inner protection sleeve (03.010.435 or 03.010.436) must be removed. The nail is then inserted through the soft and flexible outer sleeve (03.010.437S or 03.010.438S).




▲ Precaution:


The knee must remain in extension once the handle assembly has been inserted.

■ Note:

The reamer can catch on the edge of the inner protection sleeve because of the sharp angle entering the tibia (1). Stop reamer rotation and lift up on the drill to center the reamer shaft in the sleeve (2). Move the reamer in and out until the reamer comes out.

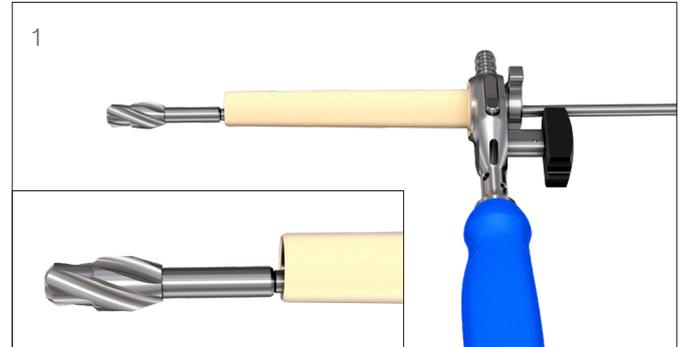

Reamer shaft angled through sleeve

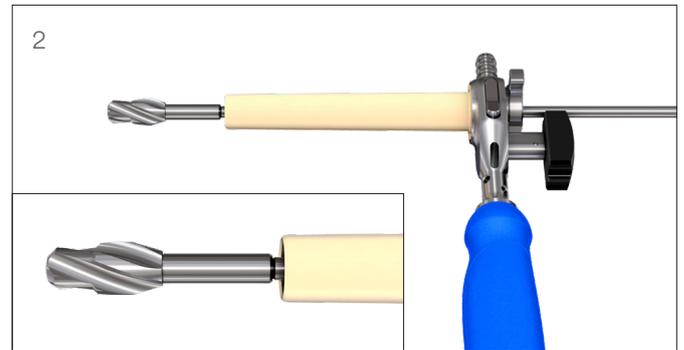

Reamer shaft centered through sleeve

Implants

Tibial Nail Advanced

Length (mm)	Ø 8 mm	Ø 9 mm	Ø 10 mm
255	04.043.005S	04.043.105S	04.043.205S
270	04.043.010S	04.043.110S	04.043.210S
285	04.043.015S	04.043.115S	04.043.215S
300	04.043.020S	04.043.120S	04.043.220S
315	04.043.025S	04.043.125S	04.043.225S
330	04.043.030S	04.043.130S	04.043.230S
345	04.043.035S	04.043.135S	04.043.235S
360	04.043.040S	04.043.140S	04.043.240S
375	04.043.045S	04.043.145S	04.043.245S
390	04.043.050S	04.043.150S	04.043.250S
405	04.043.055S	04.043.155S	04.043.255S
420	04.043.060S	04.043.160S	04.043.260S
435	04.043.065S	04.043.165S	04.043.265S
450	04.043.070S	04.043.170S	04.043.270S
465	04.043.075S	04.043.175S	04.043.275S

Length (mm)	Ø 11 mm	Ø 12 mm	Ø 13 mm
255	04.043.305S	04.043.405S	04.043.505S
270	04.043.310S	04.043.410S	04.043.510S
285	04.043.315S	04.043.415S	04.043.515S
300	04.043.320S	04.043.420S	04.043.520S
315	04.043.325S	04.043.425S	04.043.525S
330	04.043.330S	04.043.430S	04.043.530S
345	04.043.335S	04.043.435S	04.043.535S
360	04.043.340S	04.043.440S	04.043.540S
375	04.043.345S	04.043.445S	04.043.545S
390	04.043.350S	04.043.450S	04.043.550S
405	04.043.355S	04.043.455S	04.043.555S
420	04.043.360S	04.043.460S	04.043.560S
435	04.043.365S	04.043.465S	04.043.565S
450	04.043.370S	04.043.470S	04.043.570S
465	04.043.375S	04.043.475S	04.043.575S

End Cap for Tibial Nail Advanced

Article No.	Extension (mm)
04.045.850S	0
04.045.855S	5
04.045.860S	10
04.045.865S	15

Locking Screw for Medullary Nails, Ø 5 mm*

Article No.	Length (mm)	Article No.	Length (mm)
04.045.026	26	04.045.062	62
04.045.028	28	04.045.064	64
04.045.030	30	04.045.066	66
04.045.032	32	04.045.068	68
04.045.034	34	04.045.070	70
04.045.036	36	04.045.072	72
04.045.038	38	04.045.074	74
04.045.040	40	04.045.076	76
04.045.042	42	04.045.078	78
04.045.044	44	04.045.080	80
04.045.046	46	04.045.082	82
04.045.048	48	04.045.084	84
04.045.050	50	04.045.086	86
04.045.052	52	04.045.088	88
04.045.054	54	04.045.090	90
04.045.056	56	04.045.095	95
04.045.058	58	04.045.100	100
04.045.060	60		

Locking Screw for Medullary Nails, Ø 4 mm*

Article No.	Length (mm)	Article No.	Length (mm)
04.045.218	18	04.045.250	50
04.045.220	20	04.045.252	52
04.045.222	22	04.045.254	54
04.045.224	24	04.045.256	56
04.045.226	26	04.045.258	58
04.045.228	28	04.045.260	60
04.045.230	30	04.045.262	62
04.045.232	32	04.045.264	64
04.045.234	34	04.045.266	66
04.045.236	36	04.045.268	68
04.045.238	38	04.045.270	70
04.045.240	40	04.045.272	72
04.045.242	42	04.045.274	74
04.045.244	44	04.045.276	76
04.045.246	46	04.045.278	78
04.045.248	48	04.045.280	80

Locking Screw for Medullary Nails, Low Profile, Ø 5 mm*

Article No.	Length (mm)	Article No.	Length (mm)
04.045.326	26	04.045.362	62
04.045.328	28	04.045.364	64
04.045.330	30	04.045.366	66
04.045.332	32	04.045.368	68
04.045.334	34	04.045.370	70
04.045.336	36	04.045.372	72
04.045.338	38	04.045.374	74
04.045.340	40	04.045.376	76
04.045.342	42	04.045.378	78
04.045.344	44	04.045.380	80
04.045.346	46	04.045.382	82
04.045.348	48	04.045.384	84
04.045.350	50	04.045.386	86
04.045.352	52	04.045.388	88
04.045.354	54	04.045.390	90
04.045.356	56	04.045.395	95
04.045.358	58	04.045.400	100
04.045.360	60		

Locking Screw for Medullary Nails, Low Profile, Ø 4 mm*

Article No.	Length (mm)	Article No.	Length (mm)
04.045.518	18	04.045.550	50
04.045.520	20	04.045.552	52
04.045.522	22	04.045.554	54
04.045.524	24	04.045.556	56
04.045.526	26	04.045.558	58
04.045.528	28	04.045.560	60
04.045.530	30	04.045.562	62
04.045.532	32	04.045.564	64
04.045.534	34	04.045.566	66
04.045.536	36	04.045.568	68
04.045.538	38	04.045.570	70
04.045.540	40	04.045.572	72
04.045.542	42	04.045.574	74
04.045.544	44	04.045.576	76
04.045.546	46	04.045.578	78
04.045.548	48	04.045.580	80

Alternatively, the Tibial Nail Advanced implants can be applied using associated instrumentation and a set of the following compatible screw implants:

Locking Screw StarDrive Ø 5 mm (light green)*

Article No.	Length (mm)	Article No.	Length (mm)
04.005.516	26	04.005.548	58
04.005.518	28	04.005.550	60
04.005.520	30	04.005.552	62
04.005.522	32	04.005.554	64
04.005.524	34	04.005.556	66
04.005.526	36	04.005.558	68
04.005.528	38	04.005.560	70
04.005.530	40	04.005.562	72
04.005.532	42	04.005.564	74
04.005.534	44	04.005.566	76
04.005.536	46	04.005.568	78
04.005.538	48	04.005.570	80
04.005.540	50	04.005.575	85
04.005.542	52	04.005.580	90
04.005.544	54	04.005.585	95
04.005.546	56	04.005.590	100

Locking Screw for Medullary Nails, Ø 4 mm

Materials

Device(s)	Material(s)
Nails	Ti-6Al-4V ELI Polyetheretherketone (PEEK)
End Caps	Ti-6Al-7Nb (TAN) Titanium Alloy
Screws	Ti-6Al-7Nb (TAN) Titanium Alloy

Locking Screw StarDrive Ø 4 mm (dark blue)*

Article No.	Length (mm)	Article No.	Length (mm)
04.005.408	18	04.005.440	50
04.005.410	20	04.005.442	52
04.005.412	22	04.005.444	54
04.005.414	24	04.005.446	56
04.005.416	26	04.005.448	58
04.005.418	28	04.005.450	60
04.005.420	30	04.005.452	62
04.005.422	32	04.005.454	64
04.005.424	34	04.005.456	66
04.005.426	36	04.005.458	68
04.005.428	38	04.005.460	70
04.005.430	40	04.005.462	72
04.005.432	42	04.005.464	74
04.005.434	44	04.005.466	76
04.005.436	46	04.005.468	78
04.005.438	48	04.005.470	80

*All screws available non-sterile or in sterile packaging. Corresponding sterile article number with suffix "TS" for tube packaging or "S" for standard pouch packaging.

Screw length designations are defined to reflect the readings on the length measurement tools and do not necessarily correspond to the actual total length of the screw.

Alternative Instruments

Alternative Instruments	
03.010.104	03.010.101
	


Alternative Instruments	
03.045.022 for screw lengths up to 120 mm	03.010.061 for screw lengths up to 100 mm
	

Alternative Instruments	
03.043.001	393.100
	

03.045.035 + 03.045.036	351.717 + 351.719
	
	

03.043.015	03.043.016
	

03.045.018	357.399
	

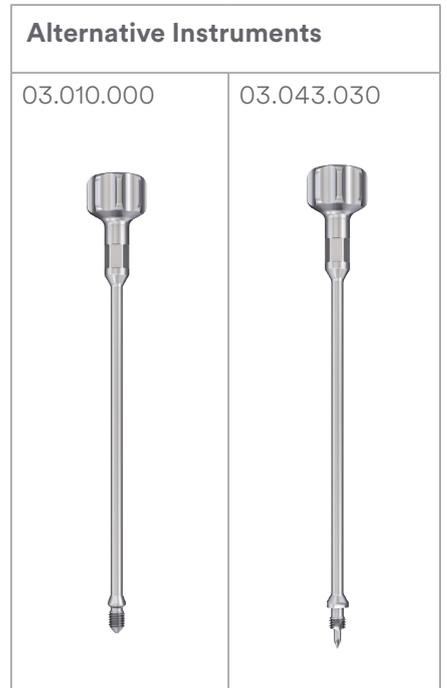

Instruments

03.010.070 Trocar Ø 4.2 mm, for No. 03.010.065



03.010.104 Drill Bit Ø 4.2 mm, calibrated, length 145 mm, 3-flute, for Quick Coupling



03.010.170 Hammer Guide



03.019.017 Depth Gauge for Locking Screws, measuring range to 100 mm



03.010.429 Direct Measuring Device for Drill Bits, length 145 mm



03.010.500 Handle, with Quick Coupling



03.010.522 Combined Hammer, 500 g



03.037.008 Awl Ø 8 / 4.7 mm, curved, cannulated



03.043.001 Universal Chuck



03.043.027 Screwdriver, with T-Handle, with spherical head, hexagonal, 8 mm



03.045.001 Screwdriver, XL25



03.045.002 Retention Pin for Screwdriver



03.045.003 Screwdriver, short, XL25



03.045.004 Retention Pin for Screwdriver, short



03.045.018	Guide Wire Ø 3.2 mm, with Drill Tip, L 400 mm	
03.045.019	Protection Sleeve, Ø 11/8	
03.045.020	Drill Sleeve, Ø 4.2 mm	
03.045.022	Drill Bit, calibrated, Ø 4.2 mm, extra-long	
03.045.035	Direct Measuring Device for Intramedullary Nails	
03.045.036	Tube for Direct Measuring Device	
321.160	Combination Wrench Ø 11.0 mm	

321.170 Pin Wrench Ø 4.5 mm, length 120 mm



03.043.003 Protection Sleeve, short



03.043.004 Wire Guide, multihole, short



03.043.016 Drill Bit, flexible, Ø 12 mm, long



03.043.028 Driving Cap



03.043.029 Aiming Arm, radiolucent



03.043.005 Wire Guide, for Nails Ø 8–11 mm, multihole, long



03.043.009 Wire Guide, for Nails Ø 8–13 mm, multihole, long



03.043.024 Insertion Handle, radiolucent, long



03.043.025 Connecting Screw, cannulated, long



03.043.026 Compression Screw, long



03.043.033S Protection Sleeve, for Nails Ø 8–11 mm, flexible, long, sterile



03.043.008S Trocar, for Nails Ø 8–11 mm, long, sterile



03.043.034S Protection Sleeve, for Nails Ø 8–13 mm, flexible, long, sterile



03.043.012S Trocar, for Nails Ø 8–13 mm, long, sterile



03.010.101 Drill Bit Ø 4.2 mm, calibrated, length 145 mm, 3-flute, with Coupling for RDL



03.010.061 Drill Bit Ø 4.2 mm, calibrated, length 340 mm, 3-flute, for Quick Coupling, for No. 03.010.065



393.100 Universal Chuck



351.717 Depth Gauge for Medullary Nails



351.719 Elongation Tube for Reaming Rods, for Depth Gauge for Medullary Nails, for Nos. 351.717 and 03.019.001



03.043.015 Drill Bit, cannulated, Ø 12 mm, long




03.043.007S Protection Sleeve, for Nails Ø 8–11 mm, flexible, Tibial Fixation, long, sterile



03.043.011S Protection Sleeve, for Nails Ø 8–13 mm, flexible, Tibial Fixation, long, sterile



03.043.035S Protection Sleeve, for Nails Ø 8–11 mm,
rigid, long, sterile

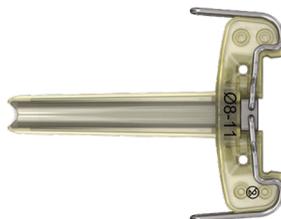

03.043.036S Protection Sleeve, for Nails Ø 8–13 mm,
rigid, long, sterile



03.043.037S Protection Sleeve, for Nails Ø 8–11 mm,
rigid, Tibial Fixation, long, sterile



03.043.038S Protection Sleeve, for Nails Ø 8–13 mm,
rigid, Tibial Fixation, long, sterile



Optional Instruments - Suprapatellar

03.010.430 Handle for Protection Sleeve



03.010.433 Centering Sleeve 12.0 / 3.2,



03.010.434 Centering Sleeve 14.5 / 3.2



03.010.435 Protection Sleeve 12.0,
for 8.0–11.0 mm Nails



03.010.436 Protection Sleeve 14.5,
for 12.0–13.0 mm Nails



03.010.437S Outer Protection Sleeve 12.0,
for 8.0–11.0 mm Nails, sterile



03.010.438S Outer Protection Sleeve 14.5,
for 12.0–13.0 mm Nails, sterile



03.010.455 Trocar Ø 12.0 mm



03.010.456 Trocar Ø 14.5 mm



Optional Instruments

03.010.093 Rod Pusher for Reaming Rod with Hexagonal Screwdriver Ø 8.0 mm



03.010.495 Intramedullary Reduction Tool, curved, with Quick Coupling, Hex 12 mm



03.010.496 T-Handle, cannulated, with Quick Coupling, Hex 12 mm



03.045.005 Screwdriver, with Quick Coupling, hexagonal 12 mm, XL25



03.045.006 Retention Pin for Screwdriver with Quick Coupling, hexagonal 12 mm



03.045.007 Screwdriver, with Quick Coupling, hexagonal 12 mm, short, XL25



03.045.008 Retention Pin for Screwdriver with Quick Coupling, hexagonal 12 mm, short



03.045.009 Sleeve for Screwdriver



03.045.010 Sleeve for Screwdriver, short



03.045.029 Reamer, Ø 5.5 mm



03.140.027 Handle, large, cannulated,
with Quick Coupling, Hex 12 mm



03.043.013 Guide Wire, with Stop, Ø 3.2 mm



Not all products are currently available in all markets.
This publication is not intended for distribution in the USA.
Please refer to the corresponding Instructions for Use for specific information on Intended use, Indications, Contraindications, Warnings and Precautions, Potential Adverse Events, Undesirable Side Effect and Residual Risks.
Instruction for Use are available at www.e-ifu.com and/or www.depuysynthes.com/ifu



Synthes GmbH
Eimattstrasse 3
4436 Oberdorf
Switzerland
Tel: +41 61 965 61 11

www.depuysynthes.com