




SPEEDTRIAD™

IMPLANT


MEDIAL TECHNIQUE OVERVIEW

1

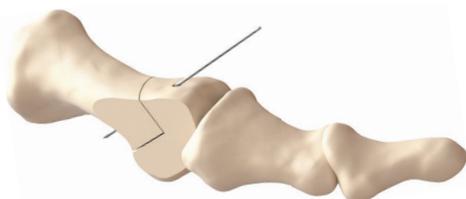


Expose the distal first metatarsal head and perform a medial exostectomy. Perform the osteotomy with the angle of the osteotomy according to surgeon preference and shift the metatarsal head laterally until the desired correction is achieved. Hold the osteotomy reduction with a temporary 0.045" K-wire and remove the proximal medial overhanging bone to create a flat surface for the implant. **NOTE:** The K-wire needs to be inserted in such a way so that it does not interfere with the fixation.

2

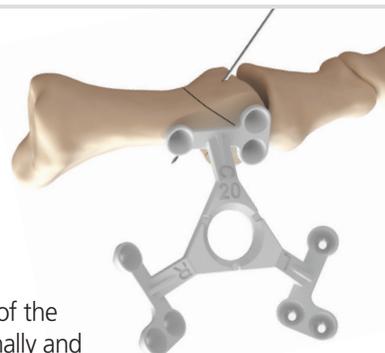


Determine the correct implant size and shape by using the BME SPEEDTRIAD™ Implant Sizing Guide (SG-2). The single leg of the implant will be placed proximally and the dual legs distally. If desired, use a sterile marker through the holes of the Sizing Guide to mark the positions of the drill holes. **NOTE:** The implant is available with centralized legs and offset distal legs.

3

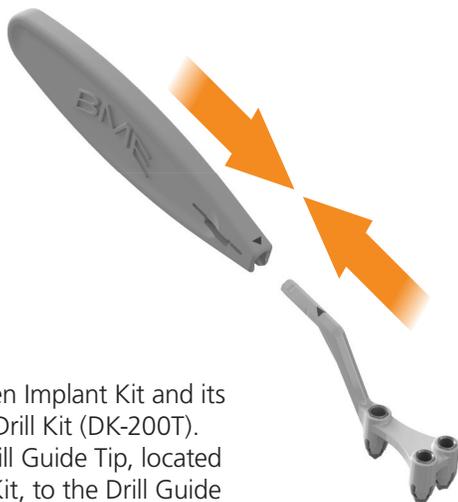


Open the chosen Implant Kit and its corresponding Drill Kit (DK-200T). Connect the Drill Guide Tip, located in the Implant Kit, to the Drill Guide Handle located in the Drill Kit.

4

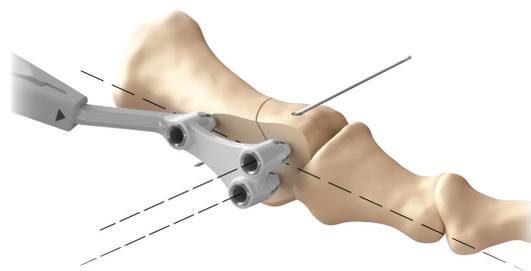


While ensuring that the bone segments are in full contact, place the Drill Guide Assembly across the osteotomy site. The implant should be positioned such that the apex of the osteotomy lies between the two distal legs. **NOTE:** The orientation of the guide should prevent inadvertent perforation of the metatarsal head into the joint when making the drill holes. The axis of the Drill Guide must be perpendicular to the osteotomy and the metatarsal head. All three prongs of the Drill Guide should be in contact with bone at all times while drilling.

5

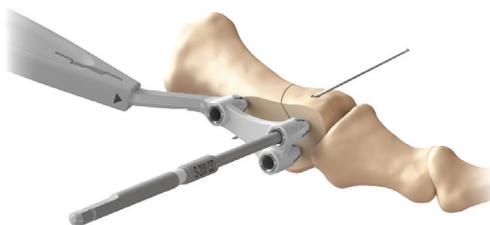

Drill the first hole on the distal side using the Drill Bit provided in the Drill Kit until the positive stop is reached.

6

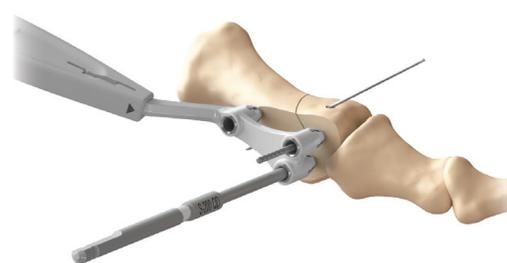


Insert a Pull Pin into the first distal hole and drill the second distal hole.

7




Insert a Pull Pin into the second distal hole, and while ensuring that the bone segments are in full contact, drill the proximal hole. The Drill Guide Assembly can be removed leaving the Pull Pins in place to mark the position of the drill holes.

8

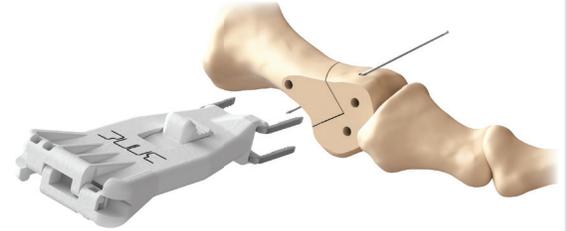


Remove the Insertion Tool containing the SPEEDTRIAD Memory Implant from the implant package. Remove the Pull Pins from the pre-drilled holes and align the tips of the legs of the SPEEDTRIAD Implant parallel with the drill holes.

9

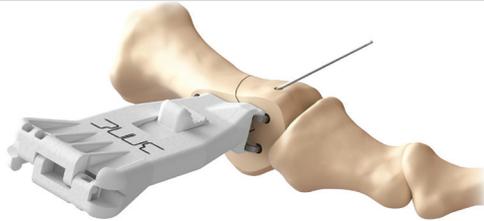


Insert the SPEEDTRIAD Memory Implant as far as possible into the pre-drilled holes. **NOTE:** To ensure proper implant placement, fluoroscopy may be used prior to releasing the implant.

10

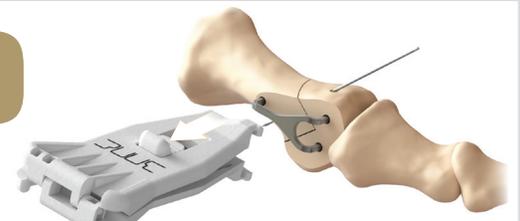


Pull the slider button away from the implant to release the implant from the Insertion Tool. **NOTE:** After pulling the slider button, move the Insertion Tool proximally to aid in release.

11

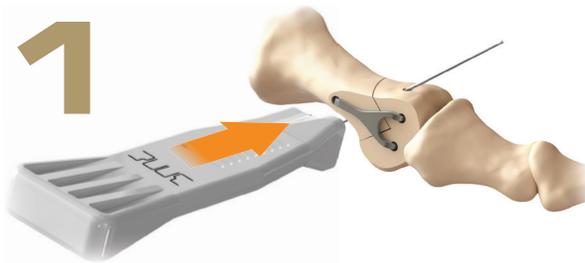

Align the supplied Tamp with the bridge of the implant, and while stabilizing the metatarsal head, use the Tamp as needed to completely seat the implant. Remove the provisional K-wire and close using established technique.

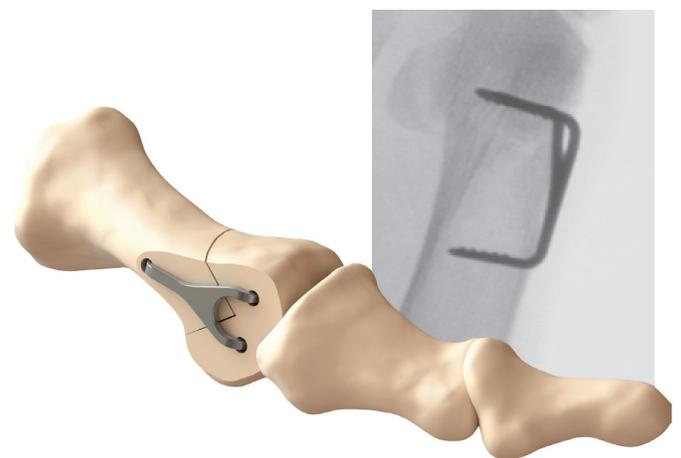

Limited Warranty and Disclaimer: DePuy Synthes Products are sold with a limited warranty to the original purchaser against defects in workmanship and materials. Any other express or implied warranties, including warranties of merchantability or fitness, are hereby disclaimed.

CAUTION: Federal Law restricts these devices to sale by or on the order of a physician. Some devices listed in this surgical technique overview may not have been licensed in accordance with Canadian law and may not be for sale in Canada.

Please contact your sales consultant for items approved for sale in Canada.

Not all Products may currently be available in all markets.

The above steps are an overview of the surgical technique. Complete information regarding indications, contra-indications, warnings, care and caution can be found in the Instructions For Use.

REMOVAL: 1. Expose the site and the bridge of the implant. 2. Using forceps, grasp the center of the implant and remove. If the implant is recessed, then use an elevator to lift the implant bridge and then use forceps to remove the implant. If solidly connected, straight implants can be removed by cutting the center of the bridge of the implant and removing the remnants with an elevator.



Manufactured or Distributed By:

Synthes USA, LLC
1101 Synthes Avenue
Monument, CO 80132

BioMedical Enterprises, Inc.
14785 Omicron Dr., Suite 205
San Antonio, TX 78245

NOTE: For recognized manufacturer, refer to the product label.

www.depuyssynthes.com

© DePuy Synthes 2016. All rights reserved.
DSUS/TRM/0716/0958(1) 10/16

TO ORDER THIS PRODUCT,
PLEASE CONTACT:



BioMedical Enterprises

(800) 880-6528 Toll Free
(866) 913-3977 Fax
(210) 677-0354 Corporate
sales@bme-tx.com Email

www.bme-tx.com